



EUROFINS PRODUCT SERVICE GMBH



TEST-REPORT

**FCC 47 CFR PART 15 SUBPART C
IC RSS 210 ISSUE 8**

**Blood Glucose Meter and Wireless Remote Controller
for Insulin Pump**

One Touch Ping Verio

FCC ID: VWT315

TEST REPORT NUMBER: G0M-1104-1064-C-1



Eurofins Product Service GmbH
Storkower Str. 38c, 15526 Reichenwalde,
Germany

Phone +49-33631-888 0
Fax +49-33631-888 660

TABLE OF CONTENTS

1	General Information	3
1.1	Notes	3
1.2	Testing laboratory	4
1.3	Details of approval holder	5
1.4	Application details	5
1.5	Acronyms and abbreviations	5
1.6	Test standards	6
1.7	Test item	6
1.8	Additional information	6
2	Technical test	7
2.1	Summary of test results	7
2.2	Test environment	7
2.3	Test equipment utilized	8
2.4	Sample emission level calculation	9
2.5	Test results	10
3	Informational Transmitter parameters	11
3.1	Occupied Bandwidth	11
4	Transmitter parameters	12
4.1	Fundamental field strength emissions	12
4.2	Emission radiated outside the specified frequency band	15
5	Receiver parameters	18
5.1	Receiver spurious emissions	18
6	Power Line parameters	20
6.1	AC power line conducted emissions	20
Annex A	Photos	21
Annex B	Transmitter occupied bandwidth	26
Annex C	Fundamental field strength emissions	29
Annex D	Transmitter radiated spurious emissions	36
Annex E	Receiver radiated spurious emissions	67
Annex F	AC Power line Conducted Emissions	74

1 General Information

1.1 Notes

The results of this test report relate exclusively to the item tested as specified in chapter "Description of test item" and are not transferable to any other test items.

Eurofins Product Service GmbH is not responsible for any generalizations and conclusions drawn from this report. Any modification of the test item can lead to invalidity of test results and this test report may therefore be not applicable to the modified test item.

The test report may only be reproduced or published in full. Reproducing or publishing extracts of the report requires the prior written approval of the Eurofins Product Service GmbH.

This document is subject to the General Terms and Conditions and the Testing and Certification System of Eurofins Product Service GmbH, available on request or accessible at www.pt.eurofins.com

Operator:

08.06.2011

B. Pudell

Date	Eurofins-Lab.	Name	Signature
------	---------------	------	-----------

Technical responsibility for area of testing:

08.06.2011

T. Jahn

Date	Eurofins	Name	Signature
------	----------	------	-----------

1.2 Testing laboratory

EUROFINS PRODUCT SERVICE GMBH
Storkower Strasse 38c
D-15526 Reichenwalde b. Berlin
Germany
Telephone :+49 33631 888 00
Telefax :+49 33631 888 660

DAKKS ACCREDITED TESTING LABORATORY
DAKKS-REGISTRATION NUMBER: D-PL-12092-01-01

RECOGNIZED NOTIFIED BODY EMC
REGISTRATION NUMBER: BNetzA-bS EMV-07/61

RECOGNIZED NOTIFIED BODY R&TTE
REGISTRATION NUMBER: BNetzA-bS-02/51-53

FCC FILED TEST LABORATORY
REG.-No. 96970

A2LA ACCREDITED TESTING LABORATORY
CERTIFICATE No. 1983.01

BLUETOOTH QUALIFICATION TEST FACILITY (BQTF)
ACCREDITED BY BLUETOOTH QUALIFICATION REVIEW BOARD

INDUSTRY CANADA FILED TEST LABORATORY
REG. NO. IC 3470

Test location, where different:

Name	: ./.
Street	: ./.
Town	: ./.
Country	: ./.
Telephone	: ./.
Fax	: ./.

1 General Information

1.1 Notes

The results of this test report relate exclusively to the item tested as specified in chapter "Description of test item" and are not transferable to any other test items.

Eurofins Product Service GmbH is not responsible for any generalizations and conclusions drawn from this report. Any modification of the test item can lead to invalidity of test results and this test report may therefore be not applicable to the modified test item.

The test report may only be reproduced or published in full. Reproducing or publishing extracts of the report requires the prior written approval of the Eurofins Product Service GmbH.

This document is subject to the General Terms and Conditions and the Testing and Certification System of Eurofins Product Service GmbH, available on request or accessible at www.pt.eurofins.com

Operator:

08.06.2011

B. Pudell



Date

Eurofins-Lab.

Name

Signature

Technical responsibility for area of testing:

08.06.2011

T. Jahn



Date

Eurofins

Name

Signature

1.6 Test standards

Technical standard : ☒ **FCC 47 CFR PART 15 SUBPART C**
☒ **IC RSS 210 ISSUE 8**

1.7 Test item

Description of test item : Blood Glucose Meter and Wireless Remote Controller for Insulin Pump
Type identification : One Touch Ping Verio
Brand Name : One Touch Ping Verio
Serial number : Unspecified
Hardware version : 2A
Software version : None
Equipment type : End product

Technical data

Radio type : Transceiver
Radio technology : Unspecified
Frequency range : 903 - 922MHz
Assigned frequency band : 902 - 928MHz
Tested frequencies : F₁ 903.0MHz
Tested frequencies : F₂ 911.5MHz
Tested frequencies : F₃ 921.3MHz
Antenna type(s) : integral
Antenna model(s) : Unspecified
Number of antennas : 1
Antenna gain(s) : Unspecified
Power supply : 3.0VDC
Duty cycle(s) : 100%
Operating mode(s) : semi duplex
Spreading technique : None
Modulation(s) : Frequency
Device classification : Portable Device (Human Body distance < 20 cm)

1.8 Additional information

None

2 Technical test

2.1 Summary of test results

No deviations from the technical specification(s) were ascertained in the course of the tests performed.



or

The deviations as specified in 2.5 were ascertained in the course of the tests performed.



2.2 Test environment

Temperature : 22 ... 26°C

Relative humidity content : 20 ... 75%

Air pressure : 86 ... 103kPa

Extreme conditions parameters:

V_{nom} : 3.0VDC
 V_{min} : N/A
 V_{max} : N/A

T_{nom} : 25°C

Other parameter: None

2.3 Test equipment utilized

Measurement Equipment List					
No.:	Measurement device:	Type:	Manufacturer:	Last Cal.	Next Cal.
ETS 0086	Semi-anechoic chamber	AC1	Frankonia	12.03.2010	12.03.2011
ETS 0271	Spectrum Analyzer	FSEK30	Rohde & Schwarz	19.03.2009	19.03.2011
ETS 0030	Biconical Antenna	HK 116	Rohde & Schwarz	10.02.2011	20.02.2012
ETS 0295	LPD Antenna	HL 223	Rohde & Schwarz	09.02.2011	09.02.2012
ETS 0018	Horn Antenna	BBHA 9120D	Schwarzbeck	26.08.2010	26.08.2011
ETS 0432	Amplifier-Matrix			02.06.2010	02.06.2012
ETS 0496	Spectrum Analyzer	FSP30	Rohde & Schwarz	26.08.2010	26.08.2011
ETS 0288	LISN	ESH2-Z5	Rohde & Schwarz	07.09.2010	07.09.2012

2.4 Sample emission level calculation

The following is a description of terms and a sample calculation, as appears in the radiated emissions data table. The numbers used in the calculation are for example only. There is no direct correlation to the specific data taken for the product described in this document:

Reading:

This is the reading obtained on the spectrum analyzer in dB μ V. Any external preamplifiers used are taken into account through internal analyzer settings.

A.F.:

This is the antenna factor for the receiving antenna. It is a conversion factor, which converts electric fields strengths to voltages, which can be measured directly on the spectrum analyzer. It is treated as a loss in dB. Cable losses have been included with the A.F. to simplify the calculations. The antenna factor is used in calculations as follows:

$$\text{Reading on Analyzer (dB}\mu\text{V)} + \text{A.F. (dB)} = \text{Net field strength (dB}\mu\text{V/m)}$$

Net:

This is the net field strength measurement (as shown above).

Limit:

This is the FCC Class B radiated emission limit (in units of dB μ V/m). The FCC limits are given in units of μ V/m. The following formula is used to convert the units of μ V/m to dB μ V/m:

$$\text{Limit (dB}\mu\text{V/m)} = 20 * \log (\mu\text{V/m})$$

Margin:

This is the margin of compliance below the FCC limit. The units are given in dB. A negative margin indicates the emission was below the limit. A positive margin indicates that the emission exceeds the limit.

Example only:

Reading	+	AF	=	Net Reading	:	Net reading - FCC limit = Margin
21.5 dB μ V		+ 26 dB		47.5 dB μ V/m		47.5 dB μ V/m - 57.0 dB μ V/m = -9.5 dB

2.5 Test results


Test case	Clause	Required	Result	Remarks
INFORMATIONAL TRANSMITTER PARAMETERS				
Occupied Bandwidth	IC RSS-Gen. 4.6.1	<input checked="" type="checkbox"/>		
TRANSMITTER PARAMETERS				
Fundamental field strength emissions	FCC § 15.249(a) IC RSS-210 A2.9(a)	<input checked="" type="checkbox"/>	PASS	
Emission radiated outside the specified frequency band	FCC § 15.249(d) FCC § 15.209 IC RSS-210 A2.9(b)	<input checked="" type="checkbox"/>	PASS	
RECEIVER PARAMETERS				
Radiated spurious emissions	IC RSS-Gen § 4.10 IC RSS-Gen § 6.1	<input checked="" type="checkbox"/>	PASS	
POWER LINE PARAMETERS				
AC power line conducted emissions	FCC § 15.207 IC RSS-Gen. 7.2.4	<input checked="" type="checkbox"/>	PASS	

3 Informational Transmitter parameters

3.1 Occupied Bandwidth

According RSS-Gen Section 4.6.1 the 99% emission bandwidth occupied by the modulated transmitted signal has to be reported as calculated or measured.

3.1.1 Measurement procedure

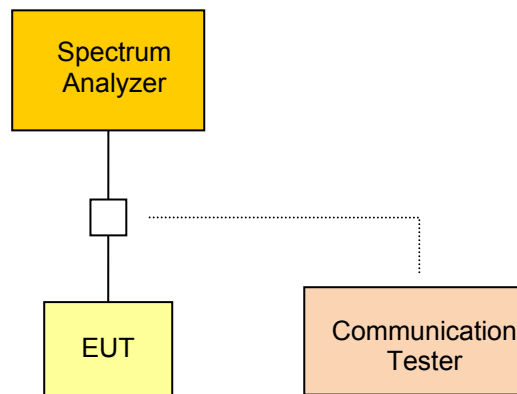

The EUT is connected to a spectrum analyzer and set to transmission mode (using a communication tester if needed) with maximum power under normal test conditions. The span of the analyzer is set wide enough to capture all significant emissions of the modulation spectrum. The resolutions bandwidth is set as close as possible to 1% of the selected span without being below 1%. The occupied bandwidth is then measured evaluated by an internal measurement procedure of the analyzer.

3.1.2 Results

Transmitter occupied bandwidth			
Measurement Conditions			
Power occupation		99%	
Channel [MHz]	Lower edge frequency [MHz]	Upper edge frequency [MHz]	Occupied Bandwidth [MHz]
903.0	902.885	903.165	0.280
911.5	911.415	911.715	0.301
921.3	921.175	921.455	0.280
See attached diagram in Annex			
Verdict			PASS

4 Transmitter parameters

4.1 Fundamental field strength emissions

According FCC rules 47 CFR 15.249(a)(1) and RSS-210 Section A2.9(a) the maximum emitted field strength at 3 meters is limited and has be verified.

4.1.1 Limits


The field strength measured at 3 metres shall not exceed the limits in the following table:

Maximum emitted field strength	
Frequency range	Fundamental field strength @ 3m
902 – 928MHz	50mV/m 94dBµV/m average 114dBµV/m peak
2400 – 2483.5MHz	50mV/m 94dBµV/m average 114dBµV/m peak
5725 – 5875MHz	50mV/m 94dBµV/m average 114dBµV/m peak

According to FCC rules the limits given in the table are average limits for frequencies above 1GHz. However, the peak field strength of any emission shall not exceed the maximum permitted average limits specified above by more than 20 dB under any condition of modulation.

According to RSS-210 the limits shown in the above table are based on measurements using an average detector, except for the fundamental emission in the frequency band 902-928 MHz, which is based on measurements using a CISPR quasi-peak detector.

4.1.2 Measurement procedure

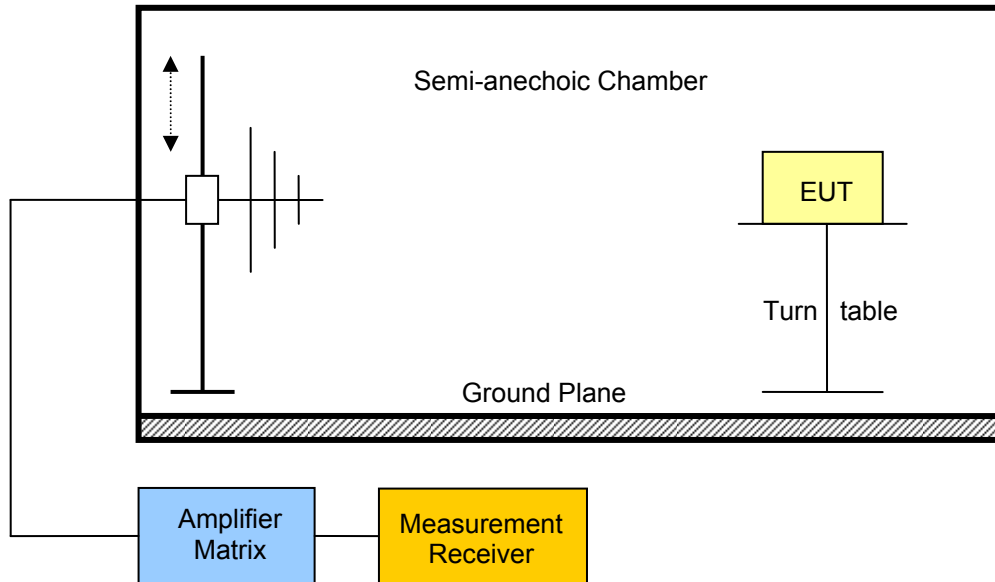

The EUT is placed on a table in a semi-anechoic chamber. The EUT is activated with the transmission modes stated in the test report. The emission level of all emission up to the 10th harmonic is scanned. In the frequency range below 1GHz a resolution bandwidth of 100kHz is used and above 1GHz a resolution bandwidth of 1MHz is used. To obtain the peak emission level the EUT is rotated through 360° and the height of the measurement antenna changed.

4.1.3 Results

Fundamental maximum field strength emissions @ 3m				
Emission [MHz]	Max. field strength [dB μ V/m]	Detector	Limit [dB μ V/m]	Margin [dB]
903.0	84.97	peak	94	-9.03
911.5	84.97	peak	94	-9.03
921.3	84.56	peak	94	-9.44
See attached diagrams in Annex				
Measurement uncertainty				4.22dB
Verdict				PASS

Comment: Due to the fact that the peak emission field-strength is below the average/quasi-peak emission limit, the corresponding average/quasi-peak measurement has been omitted and compliance with the limits is shown for the peak emissions.

The emission levels of higher order harmonics are verified in the following section.

4.2 Emission radiated outside the specified frequency band

According FCC rules 47 CFR 15.209, 15.249(d) and RSS-210 Section A2.9(b) unwanted emissions in the spurious domain are power limited and has to be validated.

4.2.1 Limits

The field strength measured at 3 metres shall not exceed the limits in the following table:

Maximum emitted field strength	
Frequency range	Field strength of harmonics
902 – 928MHz	500µV/m 54dBµV/m average 74dBµV/m peak
2400 – 2483.5MHz	500µV/m 54dBµV/m average 74dBµV/m peak
5725 – 5875MHz	500µV/m 54dBµV/m average 74dBµV/m peak

According to FCC rules the limits given in the table are average limits for frequencies above 1GHz. However, the peak field strength of any emission shall not exceed the maximum permitted average limits specified above by more than 20 dB under any condition of modulation.

According to RSS-210 the limits shown in the above table are based on measurements using an average detector, except for the fundamental emission in the frequency band 902-928 MHz, which is based on measurements using a CISPR quasi-peak detector.


Except the higher order harmonics, emission radiated outside the specified frequency band shall be attenuated by at least 50 dB below the level of the fundamental or to the general field strength limits listed in 15.209 / RSS-Gen, whichever is less stringent.

General spurious emission limits				
Frequency range [MHz]	Detector	Limit [µV/m]	Calculated Limit 3m [dBµV/m]	Measurement Distance [m]
30 – 88	Quasi-Peak	100	40	3
88 – 216	Quasi-Peak	150	43.5	3
216 – 960	Quasi-Peak	200	46	3
960 – 1000	Quasi-Peak	500	54	3
> 1000	Average	500	54	3

When average radiated emission measurements are specified, including average emission measurements below 1000 MHz, there also is a limit on the peak level of the radio frequency emissions. The limit on peak radio frequency emissions is 20 dB above the maximum permitted average emission limit applicable to the equipment under test.

Measurement procedure

The spurious emission measurement is performed on 3m a semi-anechoic test site.

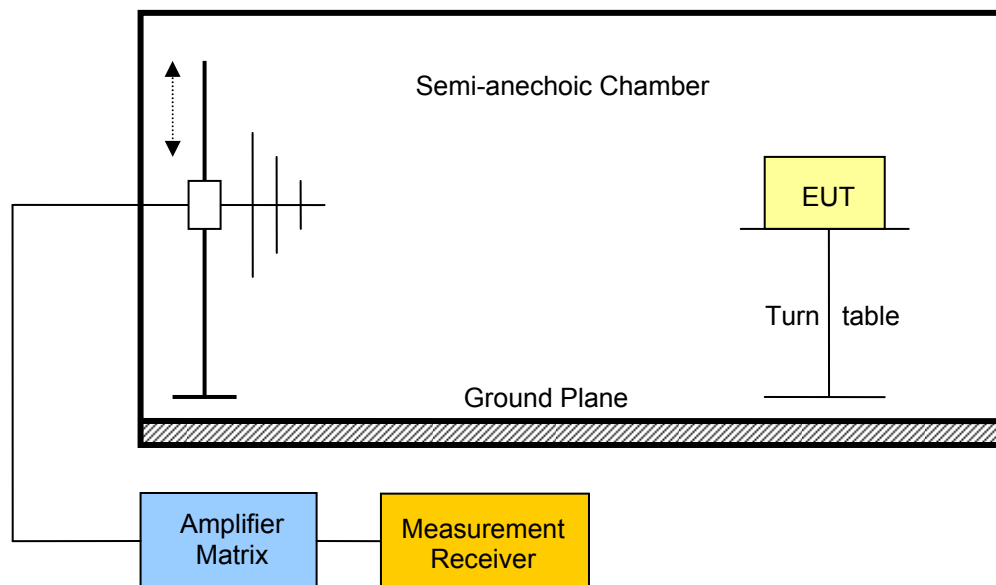

The EUT is placed on a non-metallic table. Any emission is received by the measurement antenna and measured via a measurement receiver connected to the antenna. To obtain the maximum emission the EUT is rotated through 360°.

Due to practical reasons the spurious emission level check is first performed with a peak detector and the quasi-peak and average limits.

If any emission is detected that gets close to the emission limit the detector is changed and the quasi-peak or average detector is used. Which detector is used is determined by the emission frequency. If pulsed transmission is used, averaging over the pulse train is used.

The measurement values are also corrected to obtain the field strength values at the defined measurement distances of the emission limits.

The measurement is performed over the frequency range of 30MHz up to the tenth harmonic.

4.2.2 Results

Transmitter radiated spurious emissions						
Measurement Conditions						
Measurement distance		3m				
Modulated		<input checked="" type="checkbox"/> Yes <input type="checkbox"/> No				
Channel Frequency [MHz]	Emission Frequency [MHz]	Polarization	Measured Field Strength [dBμV/m]	Limit@3m [dBμV/m]	Detector	Margin [dB]
903	902	v	41.69	46	peak	-4.31
903	902	h	44.94	46	peak	-1.06
903	5419	h	41.18	54	peak	-12.82
911	932	h	28.15	46	peak	-17.85
911	5467	h	41.44	54	peak	-12.56
921	931	v	30.89	46	peak	-15.11
921	5523	h	41.63	54	peak	-12.37
See attached diagrams in Annex						
Verdict					PASS	

Comment: Due to the fact that the peak emission field-strength is below the average/quasi-peak emission limit, the corresponding average/quasi-peak measurement has been omitted and compliance with the limits is shown for the peak emissions.

5 Receiver parameters

5.1 Receiver spurious emissions


According RSS-Gen Section 4.9 the emissions of unintentional radiators have to comply with limits stated in the rules.

5.1.1 Limits

Receiver spurious emission limits @ 3m				
Frequency range [MHz]	Detector	Limit@3m [$\mu\text{V/m}$]	Calculated Limit @ 3m [dB $\mu\text{V/m}$]	Measurement Distance [m]
30 – 88	Quasi-Peak	100	40	3
88 – 216	Quasi-Peak	150	43.5	3
216 – 960	Quasi-Peak	200	46	3
960 – 1000	Quasi-Peak	500	54	3
> 1000	Average	500	54	3

5.1.2 Measurement procedure

The spurious emission measurement is performed on a 3m test site.

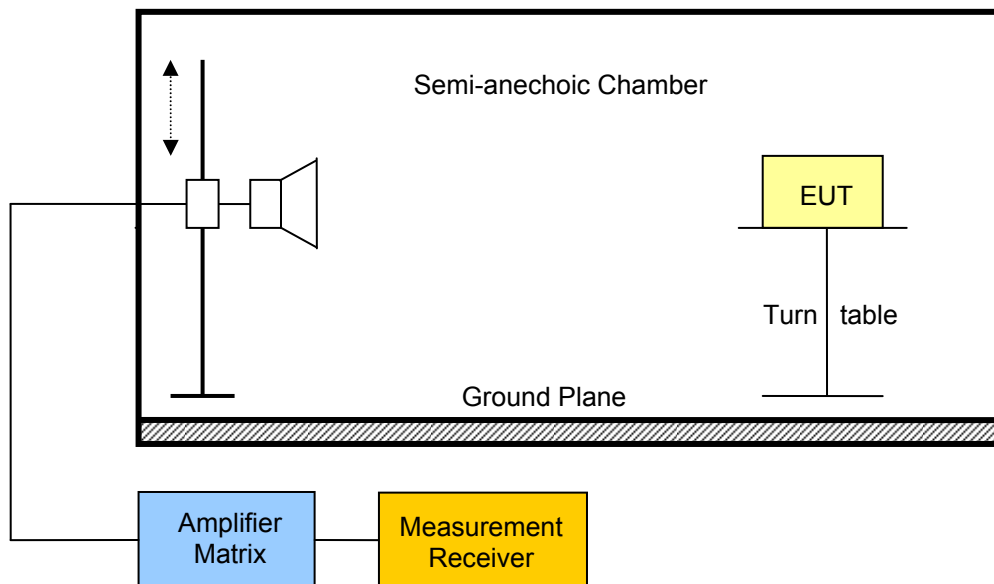

The EUT is placed on a non-metallic table. Any emission is received by a loop antenna and measured via a measurement receiver connected to the loop antenna. To obtain the maximum emission the EUT is rotated through 360°.

Due to practical reasons the spurious emission level check is first performed with a peak detector and the quasi-peak and average limits.

If any emission is detected that gets close to the emission limit the detector is changed and the quasi-peak or average detector is used. Which detector is used is determined by the emission frequency. If pulsed transmission is used, averaging over the pulse train is used.

The measurement values are also corrected to obtain the field strength values at the defined measurement distances of the emission limits.

The measurement is performed over the frequency range of 30MHz up to the 3rd harmonic.

5.1.3 Results

Receiver spurious Emissions						
Measurement Conditions						
Measurement distance		3m				
Channel Frequency [MHz]	Emission Frequency [MHz]	Polarization	Measured Field Strength [μV/m]	Limit@3m [μV/m]	Detector	Margin [μV/m]
911	897	v	50.35	200	peak	-149.65
	899	h	51.52	200	peak	-148.48
See attached diagrams in Annex						
Verdict					PASS	

Comment: Due to the fact that the peak emission field-strength is below the average/quasi-peak emission limit, the corresponding average/quasi-peak measurement has been omitted and compliance with the limits is shown for the peak emissions.

6 Power Line parameters

6.1 AC power line conducted emissions

According FCC rules 47 CFR 15.207 and RSS-Gen Section 7.2.2 for any intentional radiator that is designed to be connected to the public utility (AC) power line, the radio frequency voltage that is conducted back onto the AC power line on any frequency or frequencies within the band 150 kHz to 30 MHz shall not exceed the limits given below.

6.1.1 Limits

AC power line emission limits		
Frequency [MHz]	Conducted Limit [dB μ V]	
	Quasi-Peak	Average
0.15 – 0.5	66 to 56	56 to 46
0.5 - 5	56	46
5 - 30	60	50

6.1.2 Measurement procedure

The ac power line emissions are measured using a 50 μ H / 50 Ω line impedance stabilization network (LINS). The radio frequency voltage between each power line and ground at the power terminal is measured.


6.1.3 Results

AC power line emissions	
Conducted emission level	
See attached Diagram	
Verdict	PASS

Annex B Transmitter occupied bandwidth

RSS Gen Occupied Bandwidth

EUT Blood glucose meter
Model Titan
Approval Holder LifeScan Scotland Ltd. / Ord.: G0M-1104-1064
Temperature / Voltage tnom / Vnom
Test Site / Operator Eurofins Product Service GmbH / Mr. Treffke
Test Specification 4.4.1 Occupied Bandwidth
Comment 1 Channel.: 1
Comment 2 A spectrum analyzer with an integrated 99% power bandwidth function is used
Comment 3 OBW=280.5kHz

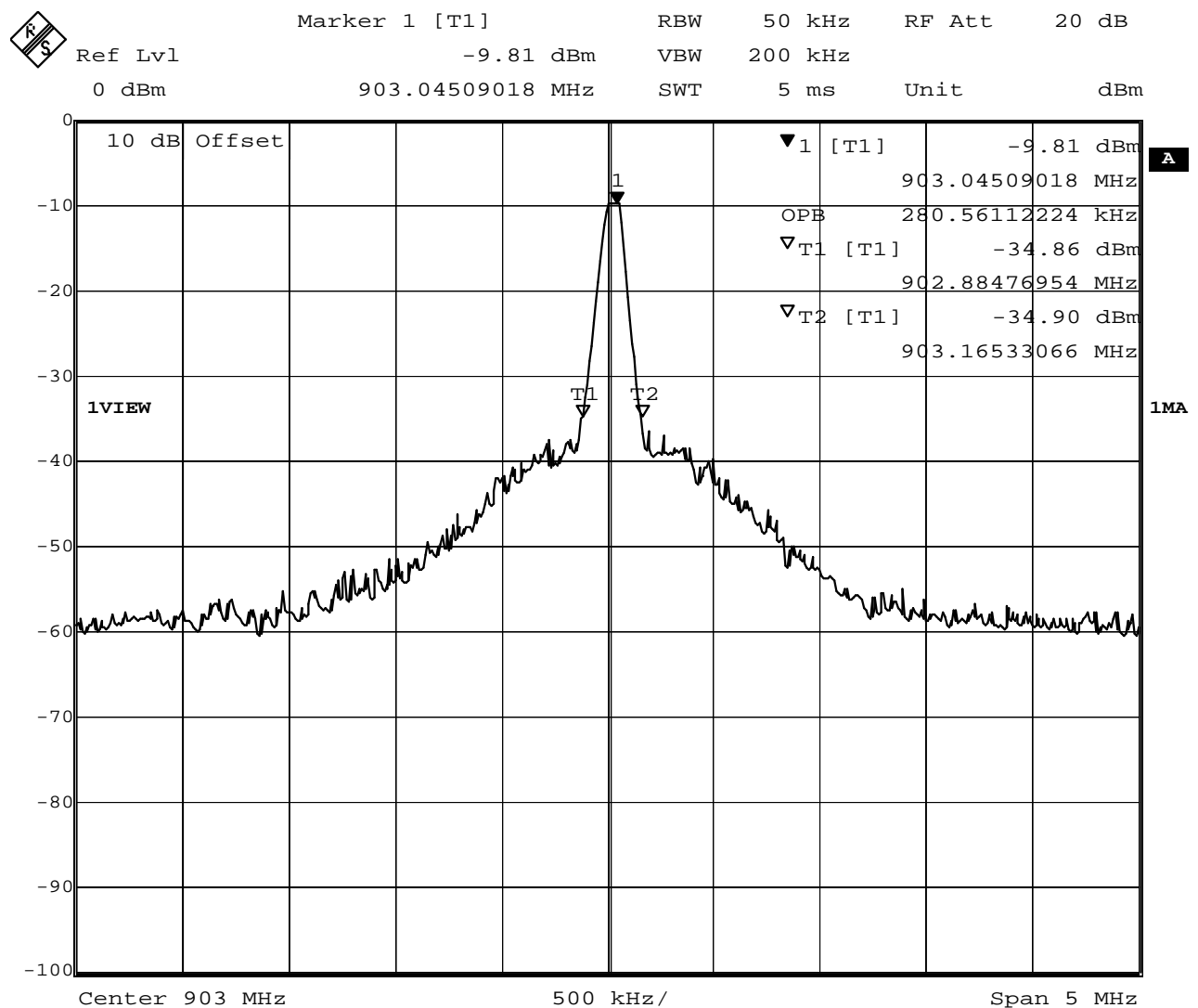

Comment A: Occupied bandwidth: 280.6 KHz

Date: 20.MAY.2011 08:47:12


Test Report No.: G0M-1104-1064-C-1

Eurofins Product Service GmbH
Storkower Str. 38c, D-15526 Reichenwalde, Germany

Page 26 of 75

RSS Gen Occupied Bandwidth

EUT Blood glucose meter
Model Titan
Approval Holder LifeScan Scotland Ltd. / Ord.: G0M-1104-1064
Temperature / Voltage tnom / Vnom
Test Site / Operator Eurofins Product Service GmbH / Mr. Treffke
Test Specification 4.4.1 Occupied Bandwidth
Comment 1 Channel.: 8
Comment 2 A spectrum analyzer with an integrated 99% power bandwidth function is used
Comment 3 OBW=300.6kHz

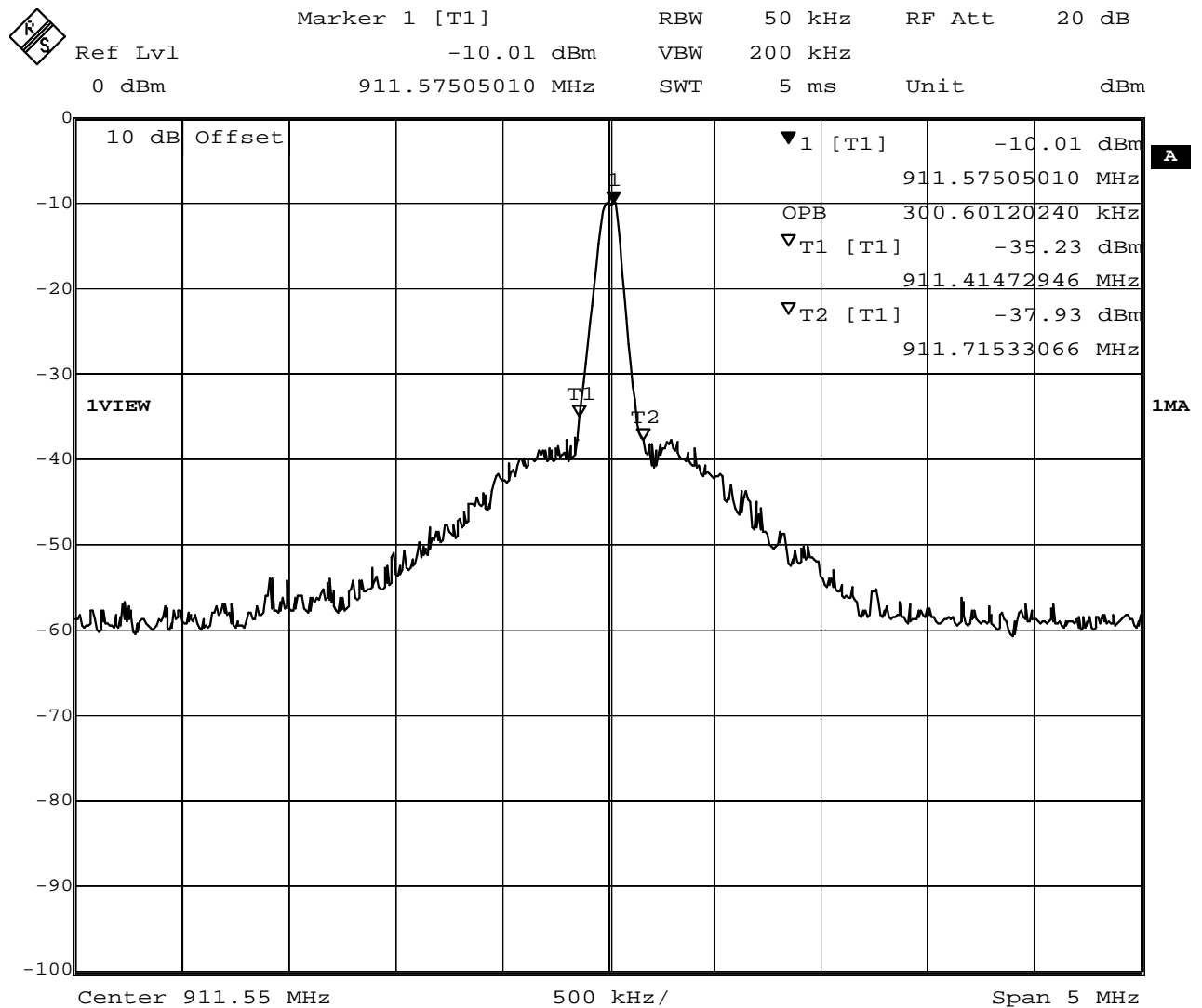

Comment A: Occupied bandwidth: 300.6 KHz
Date: 20.MAY.2011 08:53:15


Test Report No.: G0M-1104-1064-C-1

Eurofins Product Service GmbH
Storkower Str. 38c, D-15526 Reichenwalde, Germany

Page 27 of 75

RSS Gen Occupied Bandwidth

EUT Blood glucose meter
Model Titan
Approval Holder LifeScan Scotland Ltd. / Ord.: G0M-1104-1064
Temperature / Voltage tnom / Vnom
Test Site / Operator Eurofins Product Service GmbH / Mr. Treffke
Test Specification 4.4.1 Occupied Bandwidth
Comment 1 Channel.: 16
Comment 2 A spectrum analyzer with an integrated 99% power bandwidth function is used
Comment 3 OBW=280kHz

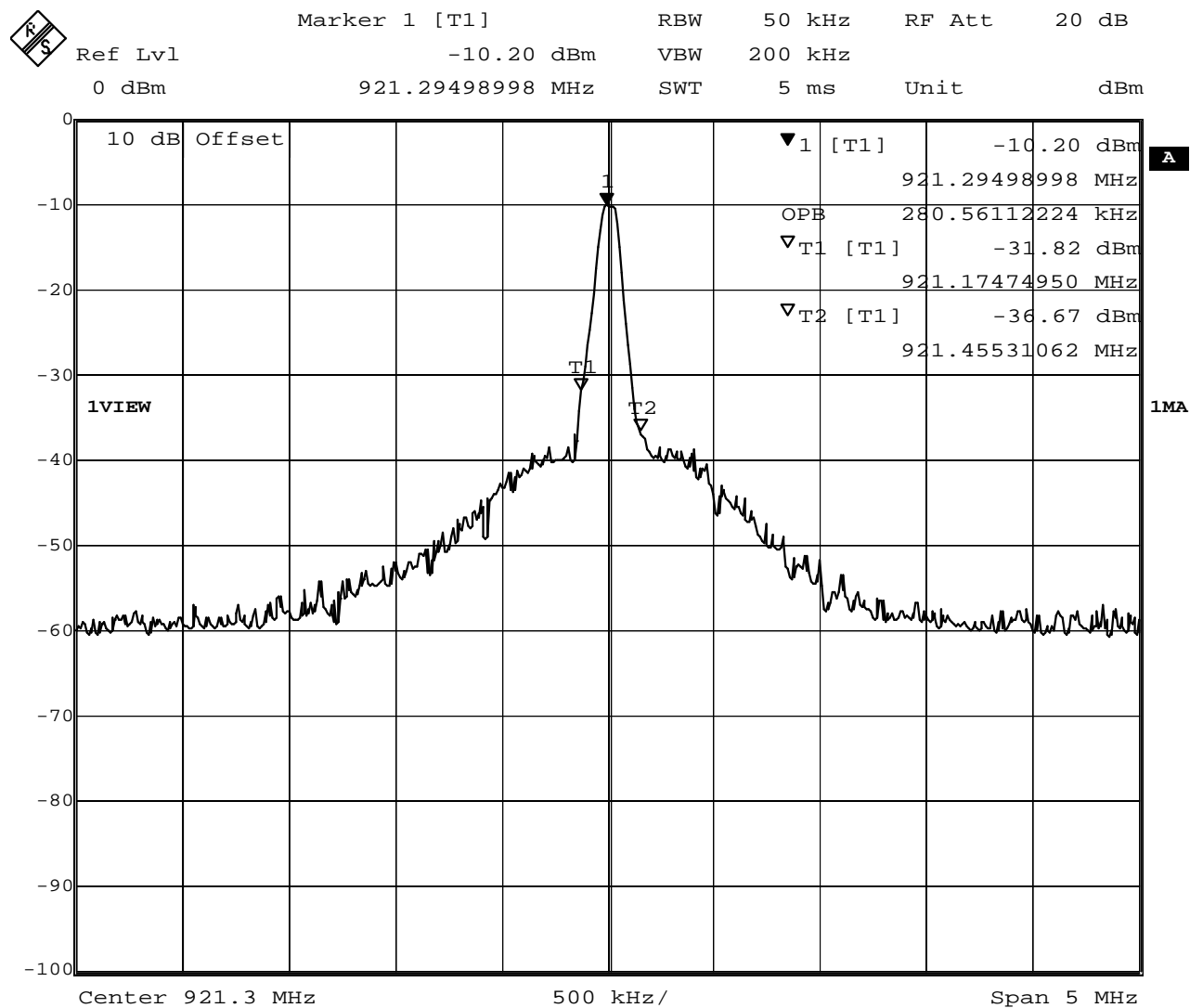


Comment A: Occupied bandwidth: 280.6 KHz
Date: 20.MAY.2011 08:51:15

Annex C Fundamental field strength emissions

Carrier power (Field Strength)

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / G0M-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 1 / Freq: 903MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 903
Comment 1: Dist.: 3m, Ant.: HL 223
Comment 2: Freq: 902.990MHz, Emax: 83.89dBµV/m, RBW: 100kHz

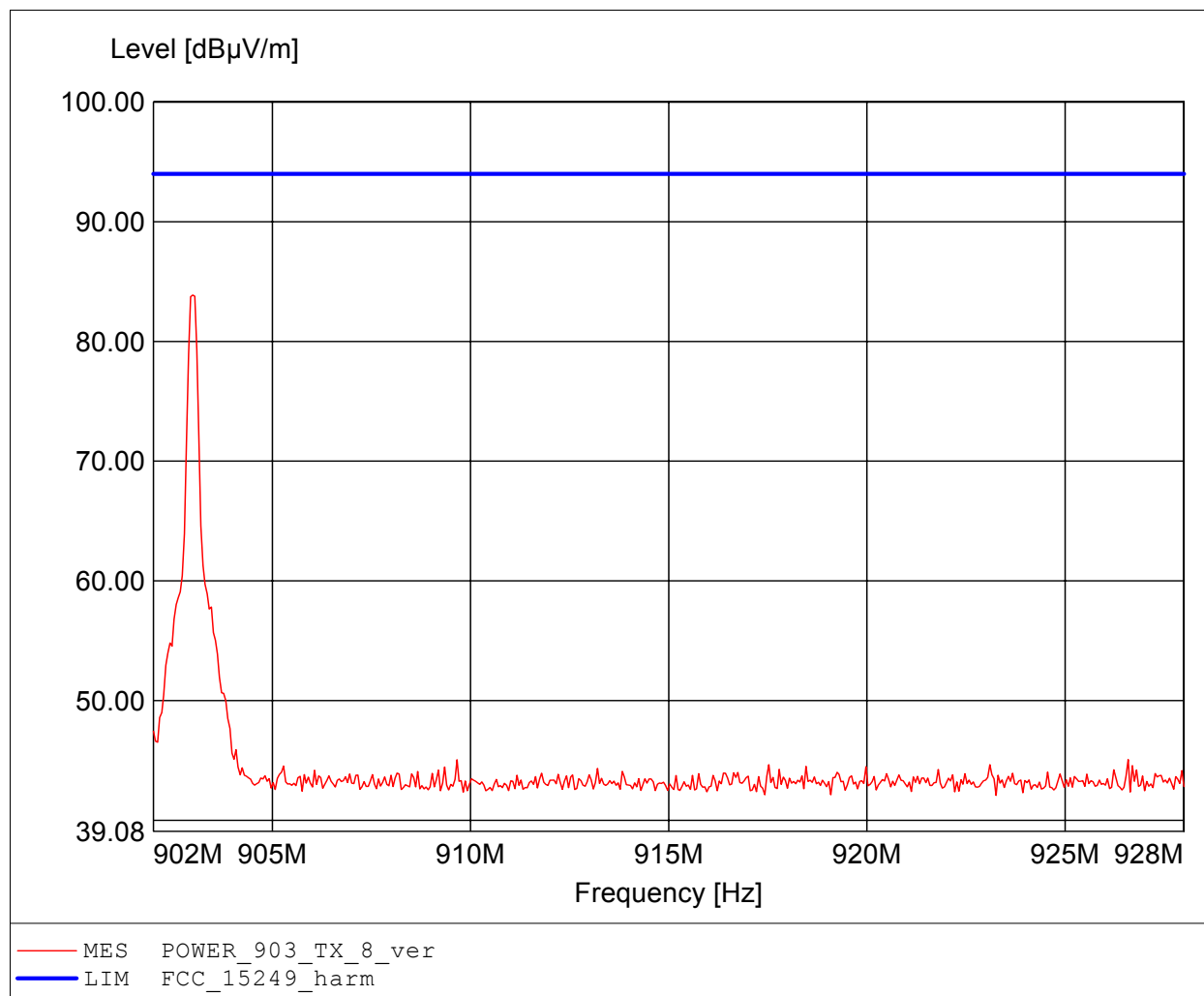

Carrier power (Field Strength)

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / G0M-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 1 / Freq: 903MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 903
Comment 1: Dist.: 3m, Ant.: HL 223
Comment 2: Freq: 902.990MHz, Emax: 84.97dBµV/m, RBW: 100kHz

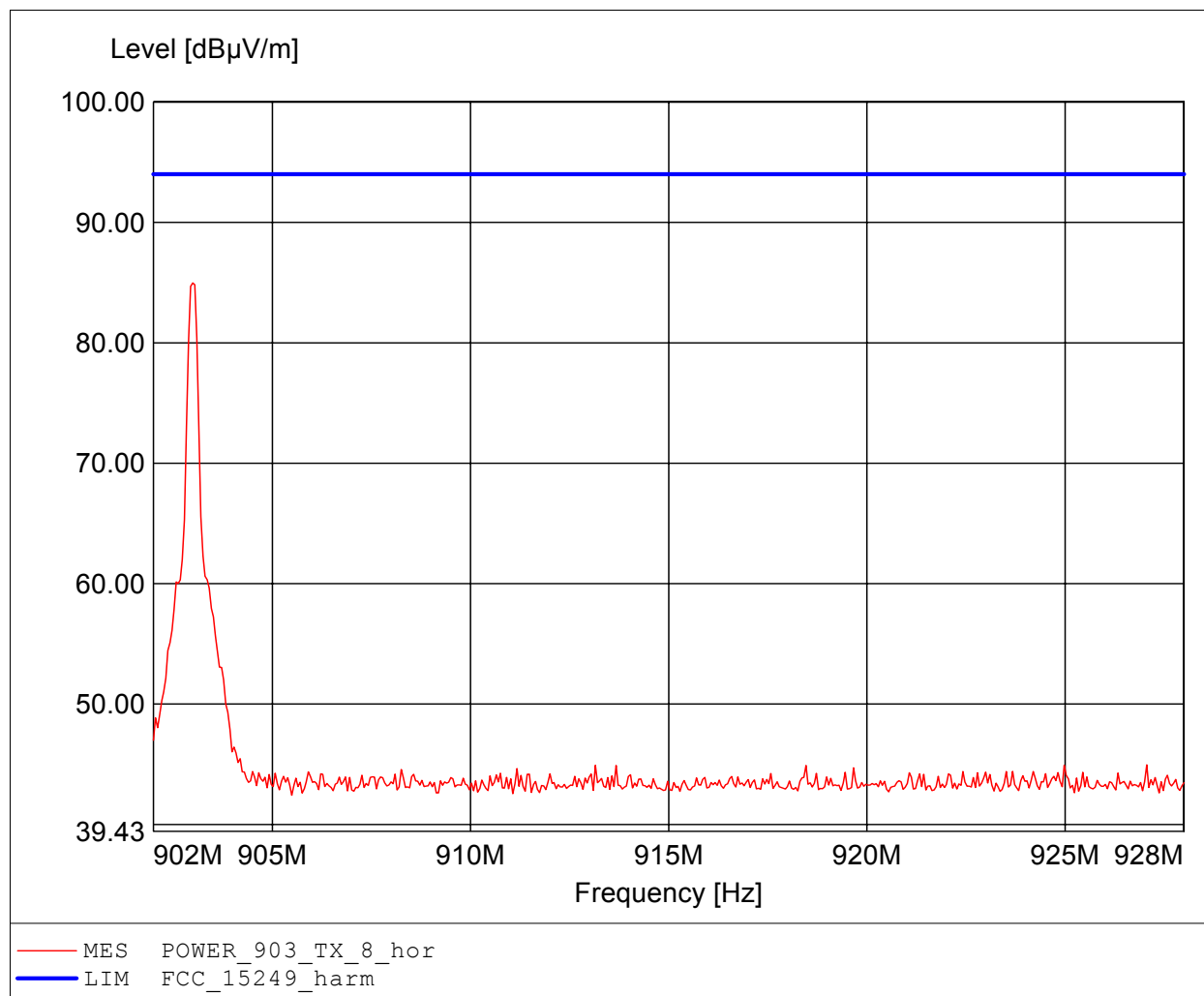

Carrier power (Field Strength)

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / GOM-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 8 / Freq: 911MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 911
Comment 1: Dist.: 3m, Ant.: HL 223
Comment 2: Freq: 911.535MHz, Emax: 84.13dB μ V/m, RBW: 1MHz

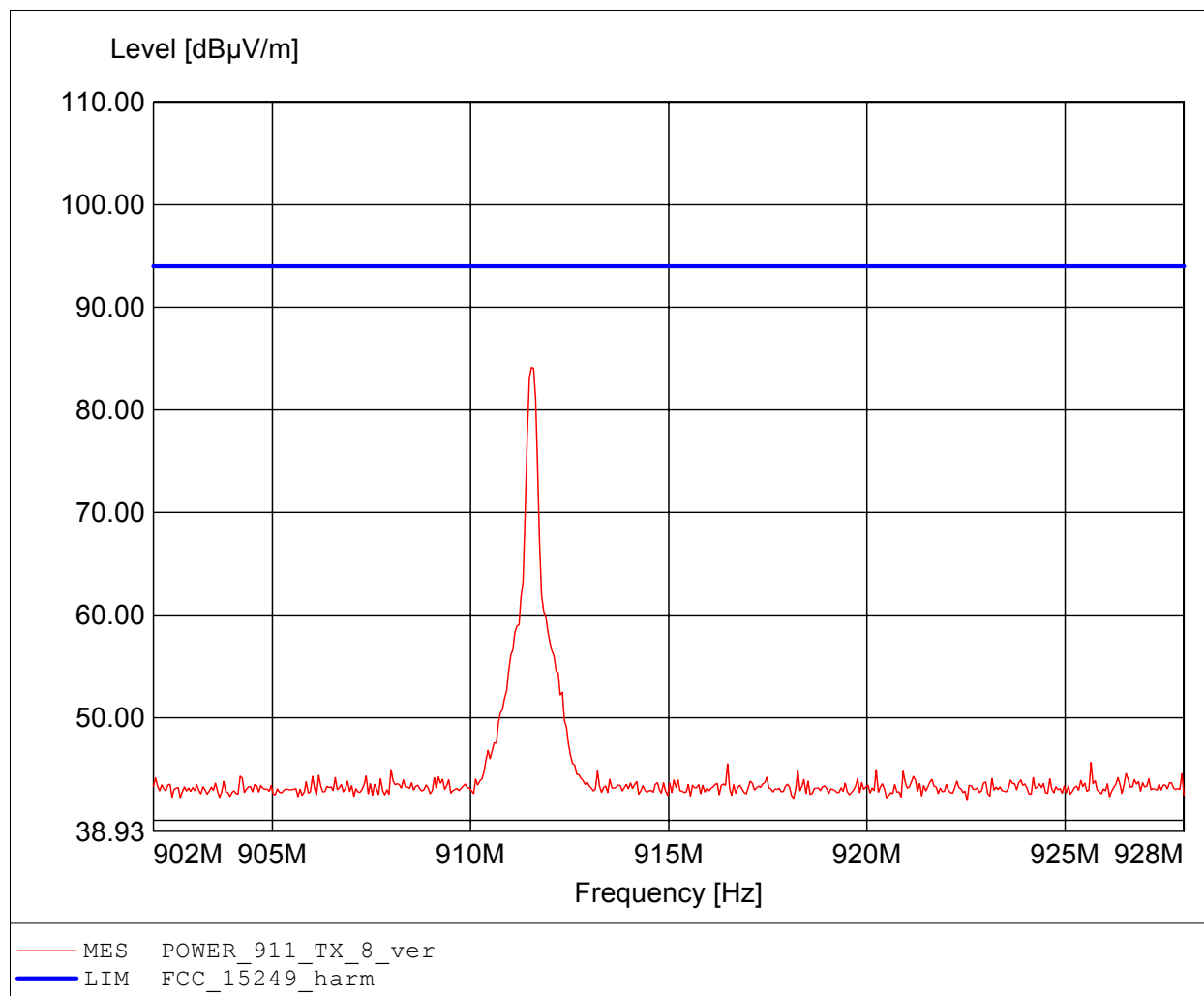

Carrier power (Field Strength)

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / G0M-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 8 / Freq: 911MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 911
Comment 1: Dist.: 3m, Ant.: HL 223
Comment 2: Freq: 911.535MHz, Emax: 84.97dBµV/m, RBW: 1MHz

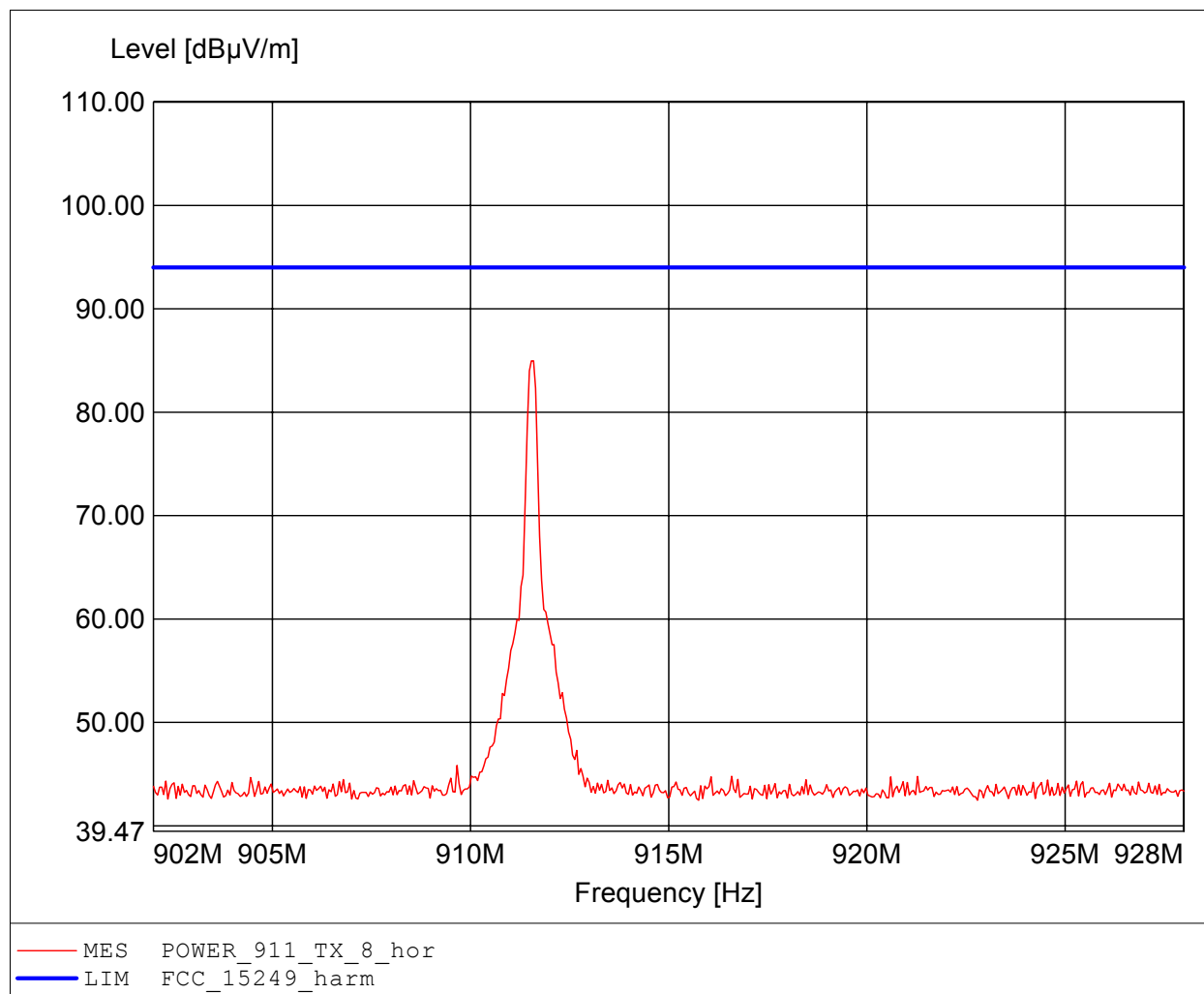

Carrier power (Field Strength)

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / GOM-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 16 / Freq: 921MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 921
Comment 1: Dist.: 3m, Ant.: HL 223
Comment 2: Freq: 921.279MHz, Emax: 83.89dBµV/m, RBW: 100kHz

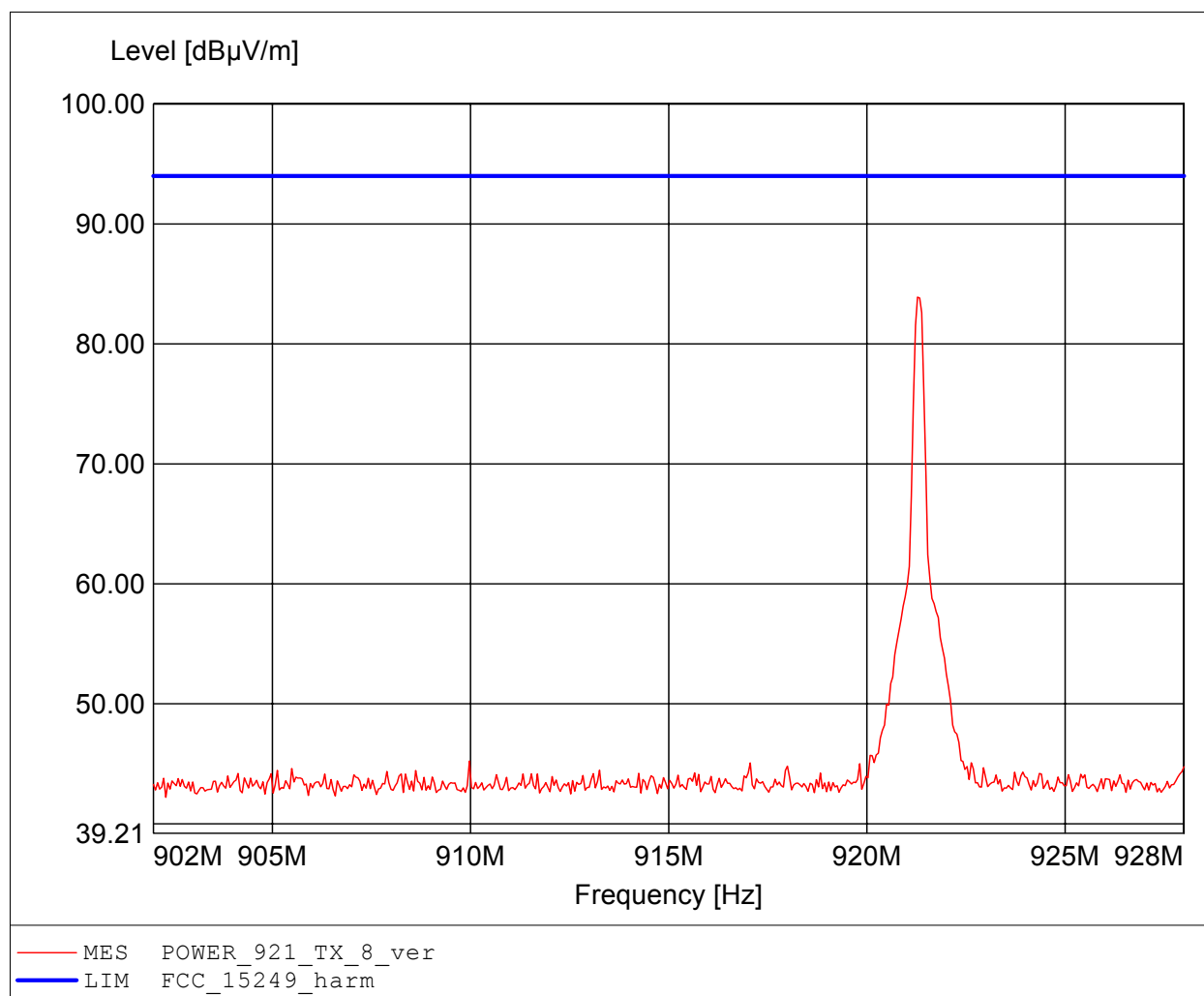

Carrier power (Field Strength)

FCC RULES PART 15, SUBPART C

Approval Holder: LifeScan Scotland Ltd. / G0M-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 16 / Freq: 921MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 921
Comment 1: Dist.: 3m, Ant.: HL 223
Comment 2: Freq: 921.279MHz, Emax: 84.56dBµV/m, RBW: 100kHz

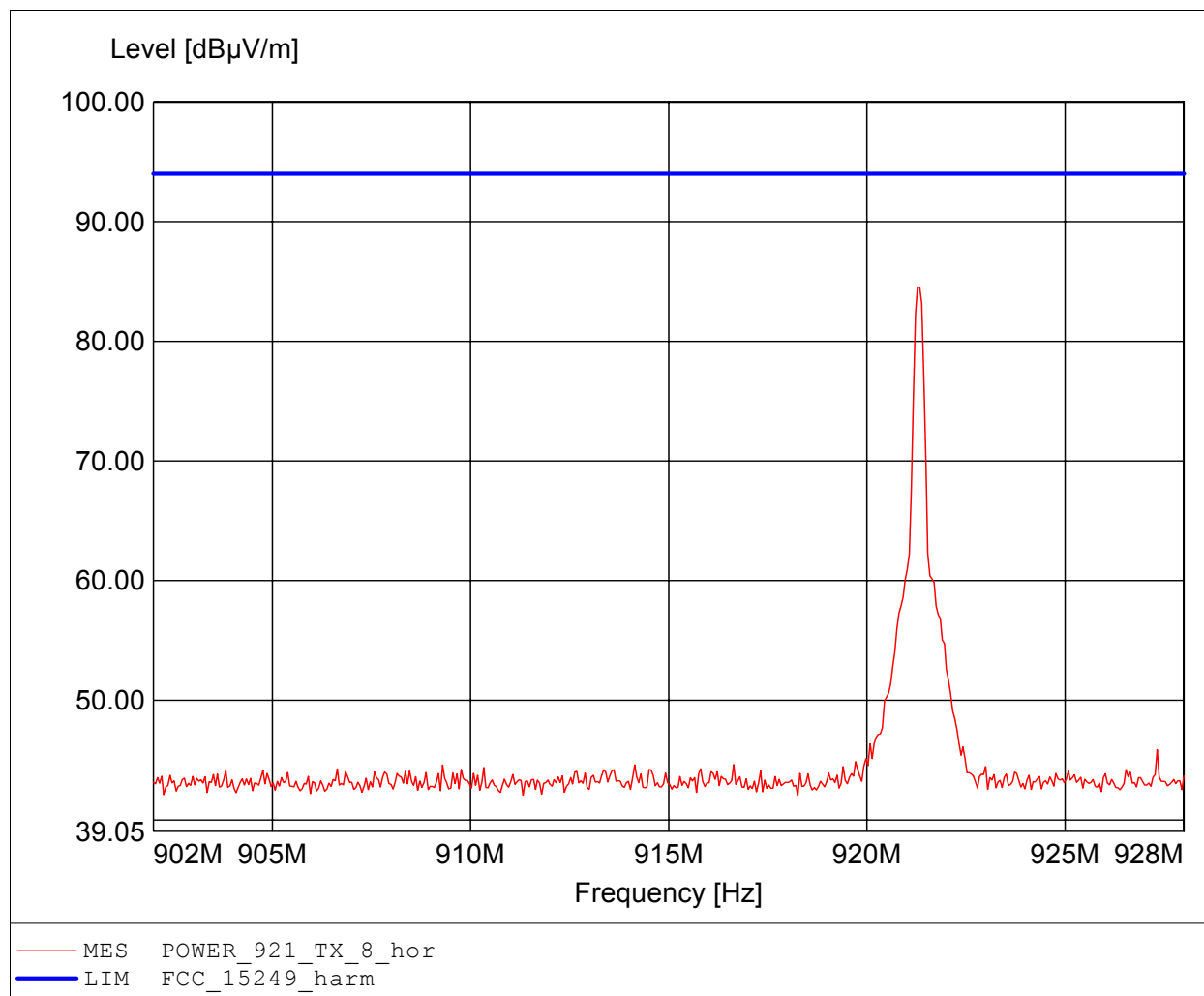


Annex D Transmitter radiated spurious emissions

Spurious emissions Field Strength

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / G0M-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 1 / Freq: 903MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 903
Comment 1: Dist.: 3m, Ant.: HK 116
Comment 2: Freq: 194.890MHz, Emax: 33.22dBµV/m, RBW: 100kHz

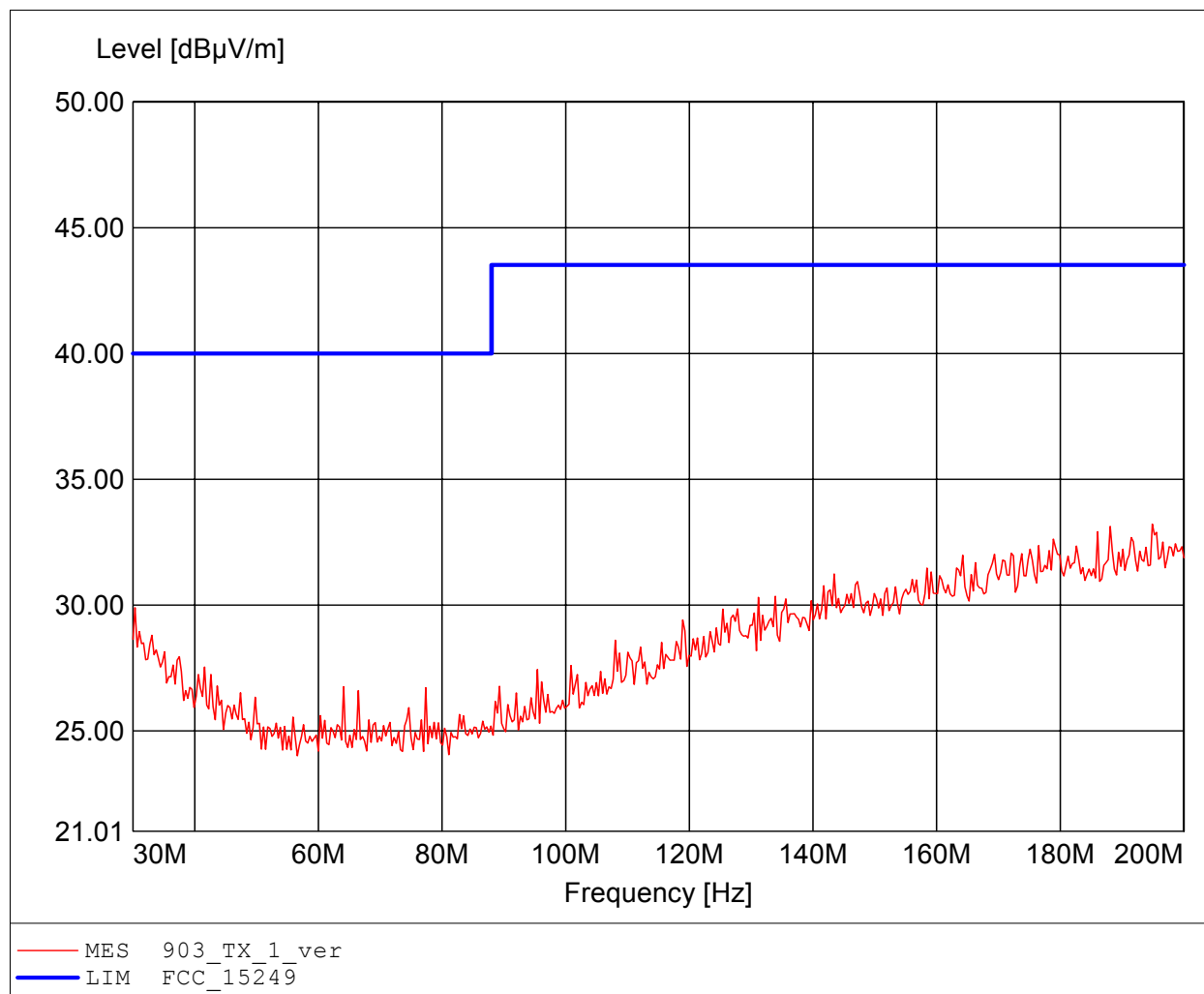

Spurious emissions Field Strength

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / G0M-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 1 / Freq: 903MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 903
Comment 1: Dist.: 3m, Ant.: HK 116
Comment 2: Freq: 192.505MHz, Emax: 33.22dBµV/m, RBW: 100kHz

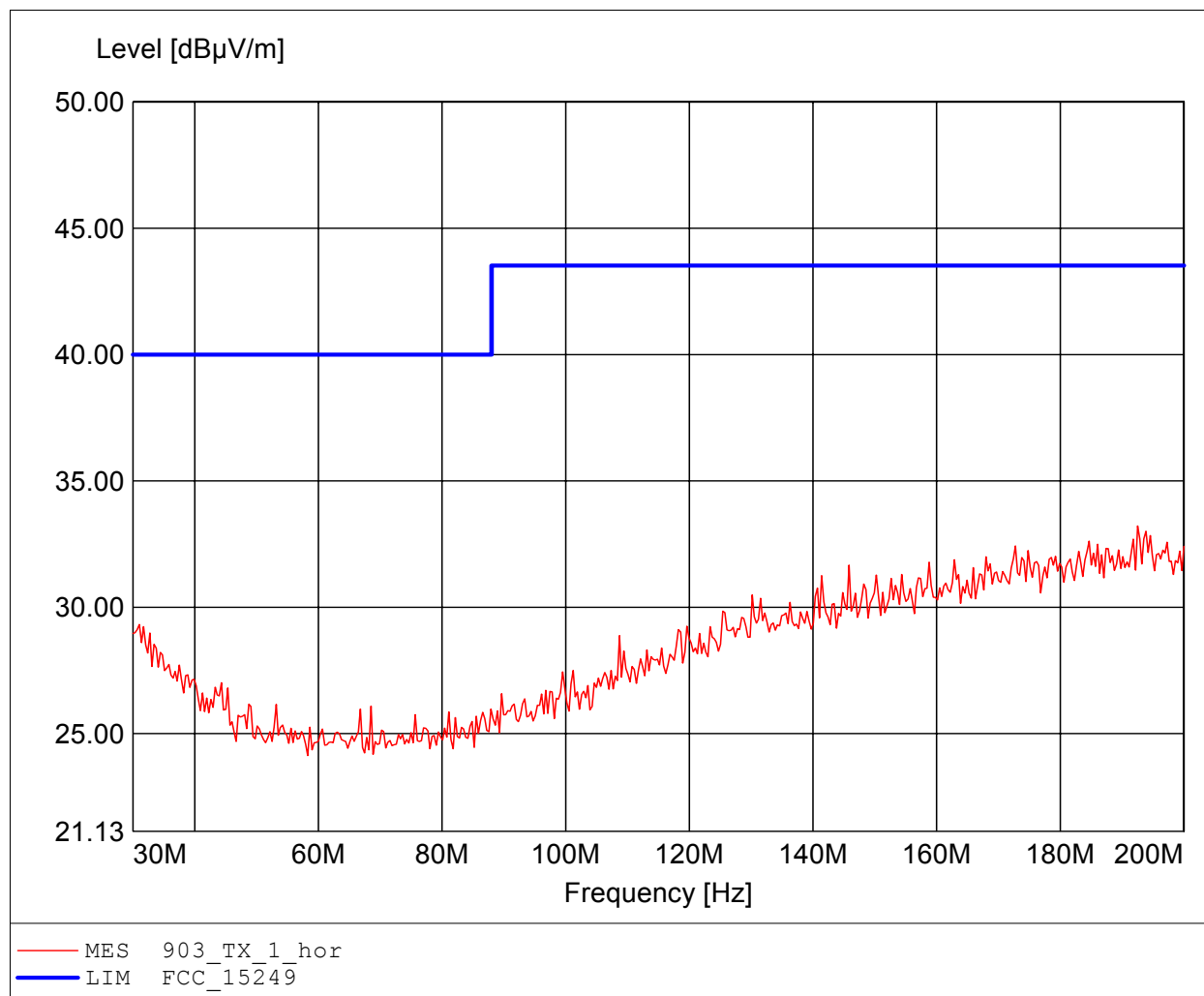

Spurious emissions Field Strength

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / G0M-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 1 / Freq: 903MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 903
Comment 1: Dist.: 3m, Ant.: HL 223, amplif.
Comment 2: Freq: 902.000MHz, Emax: 41.69dBµV/m, RBW: 100kHz

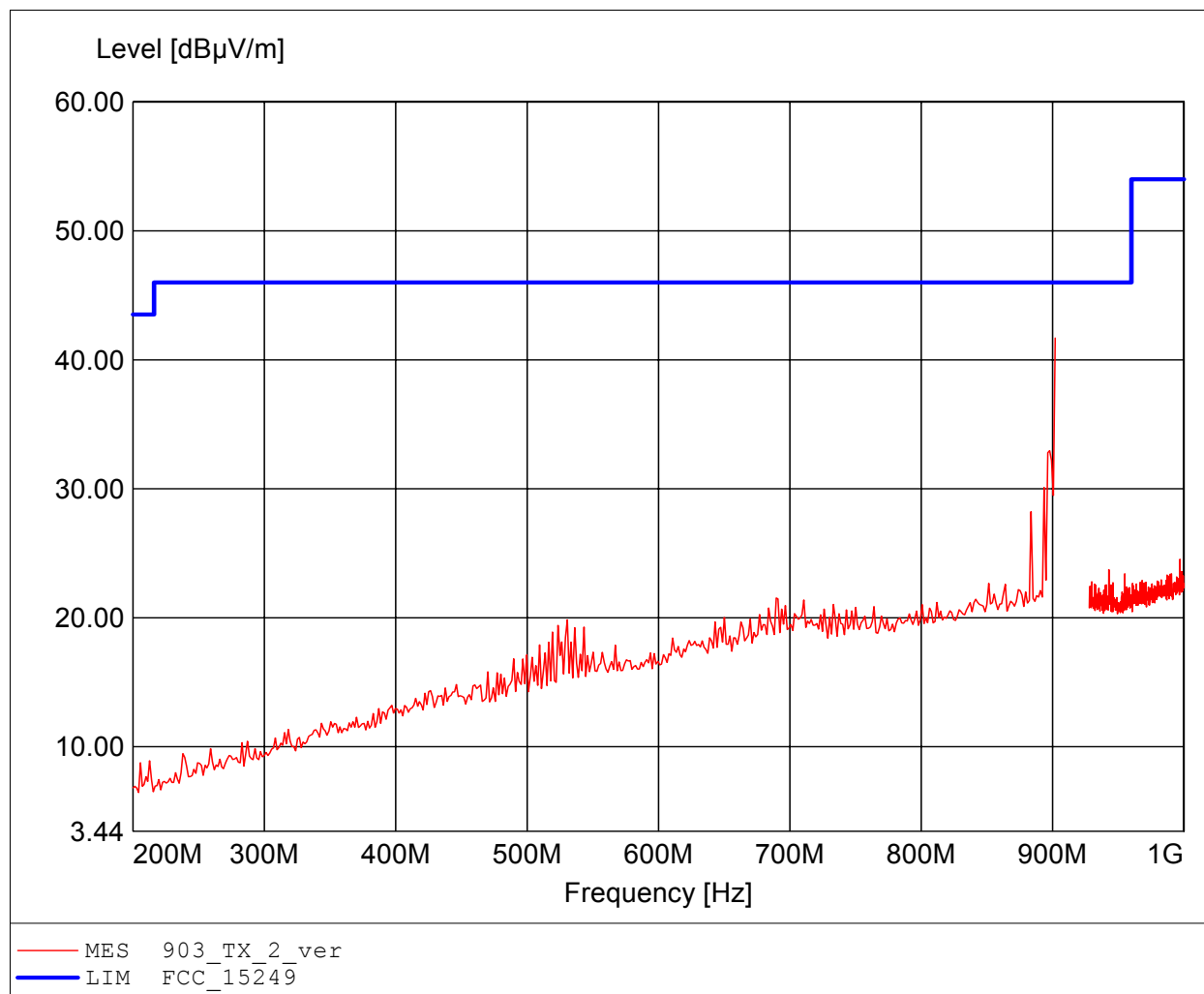

Spurious emissions Field Strength

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / G0M-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 1 / Freq: 903MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 903
Comment 1: Dist.: 3m, Ant.: HL 223, amplif.
Comment 2: Freq: 902.000MHz, Emax: 44.94dBµV/m, RBW: 100kHz

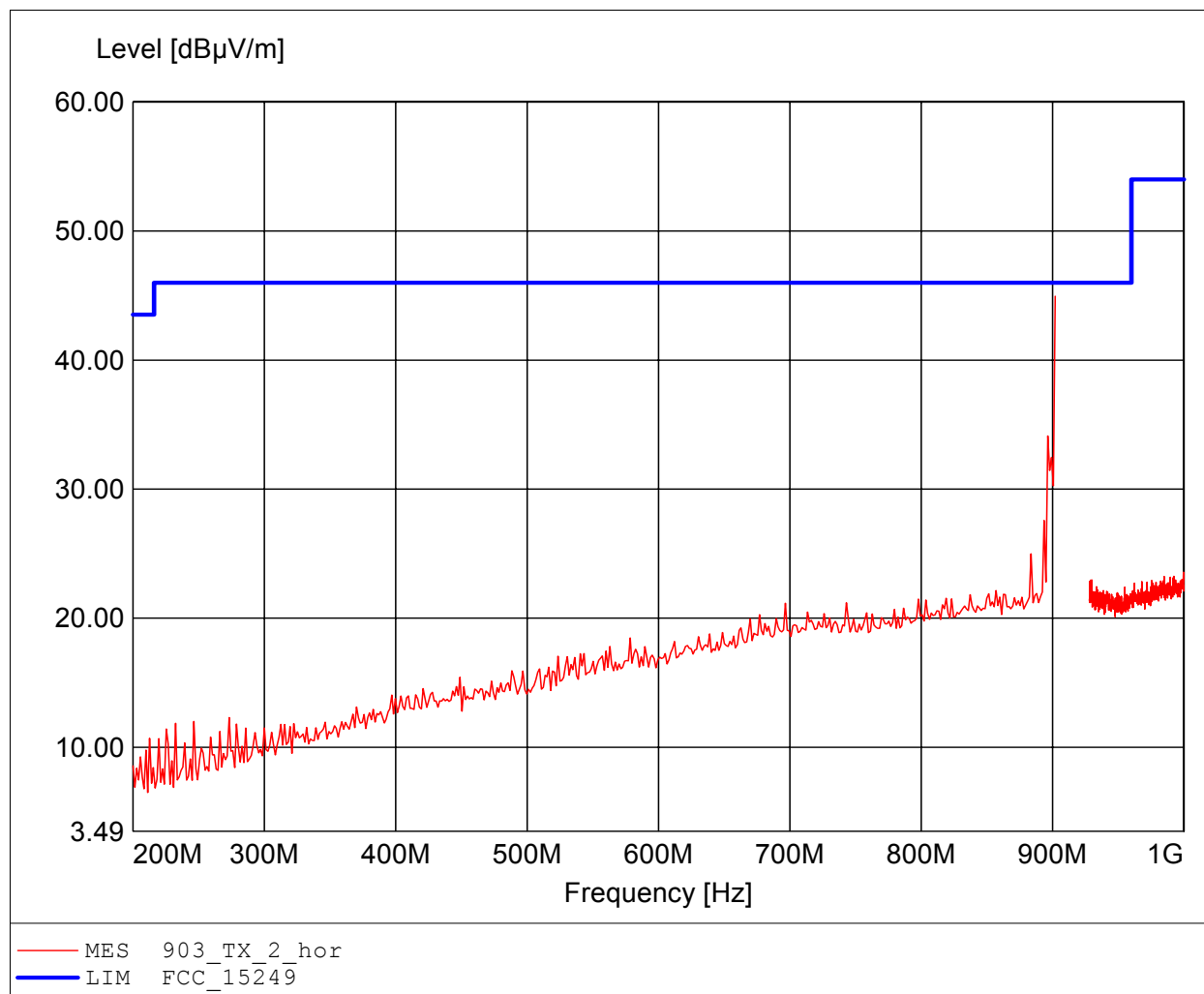

Spurious emissions Field Strength

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / GOM-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 1 / Freq: 903MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 903
Comment 1: Dist.: 3m, Ant.: BBHA9120D, amplif.
Comment 2: Freq: 3.904GHz, Emax: 42.59dB μ V/m, RBW: 1MHz

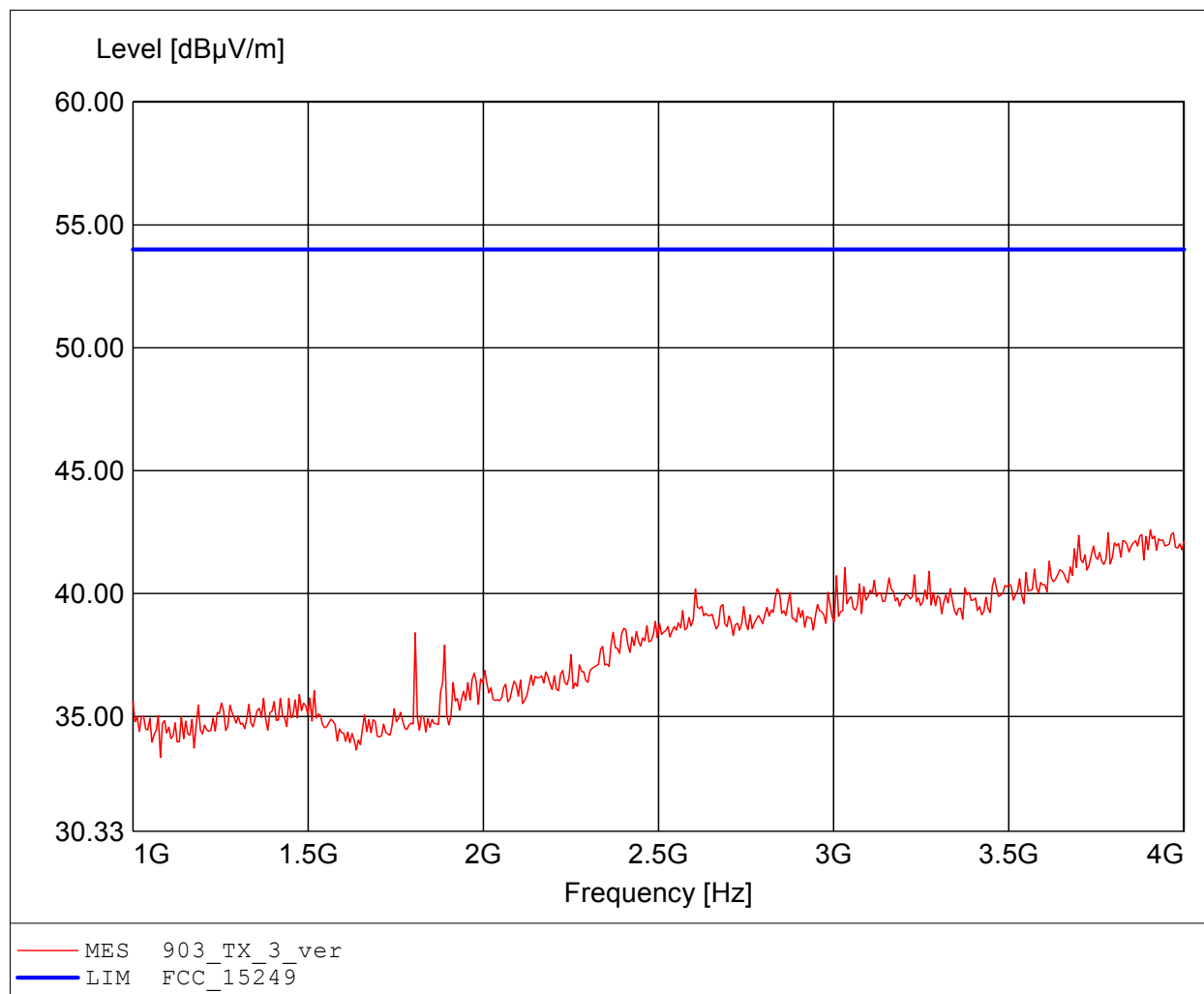

Spurious emissions Field Strength

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / G0M-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 1 / Freq: 903MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 903
Comment 1: Dist.: 3m, Ant.: BBHA9120D, amplif.
Comment 2: Freq: 3.916GHz, Emax: 43.08dBμV/m, RBW: 1MHz

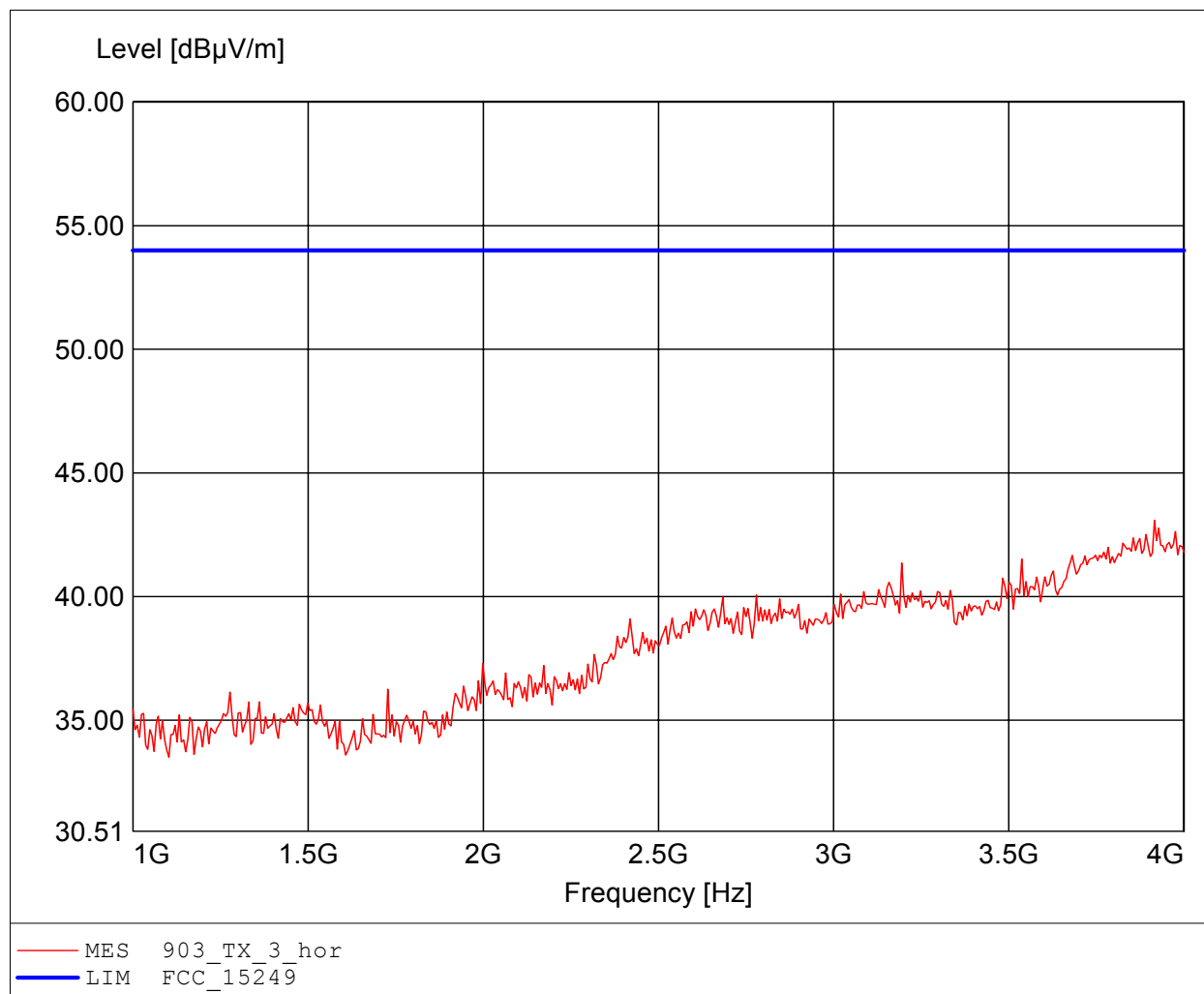

Spurious emissions Field Strength

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / GOM-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 1 / Freq: 903MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 903
Comment 1: Dist.: 3m, Ant.: BBHA9120D, ampl.+HP.
Comment 2: Freq: 7.984GHz, Emax: 40.06dBµV/m, RBW: 1MHz

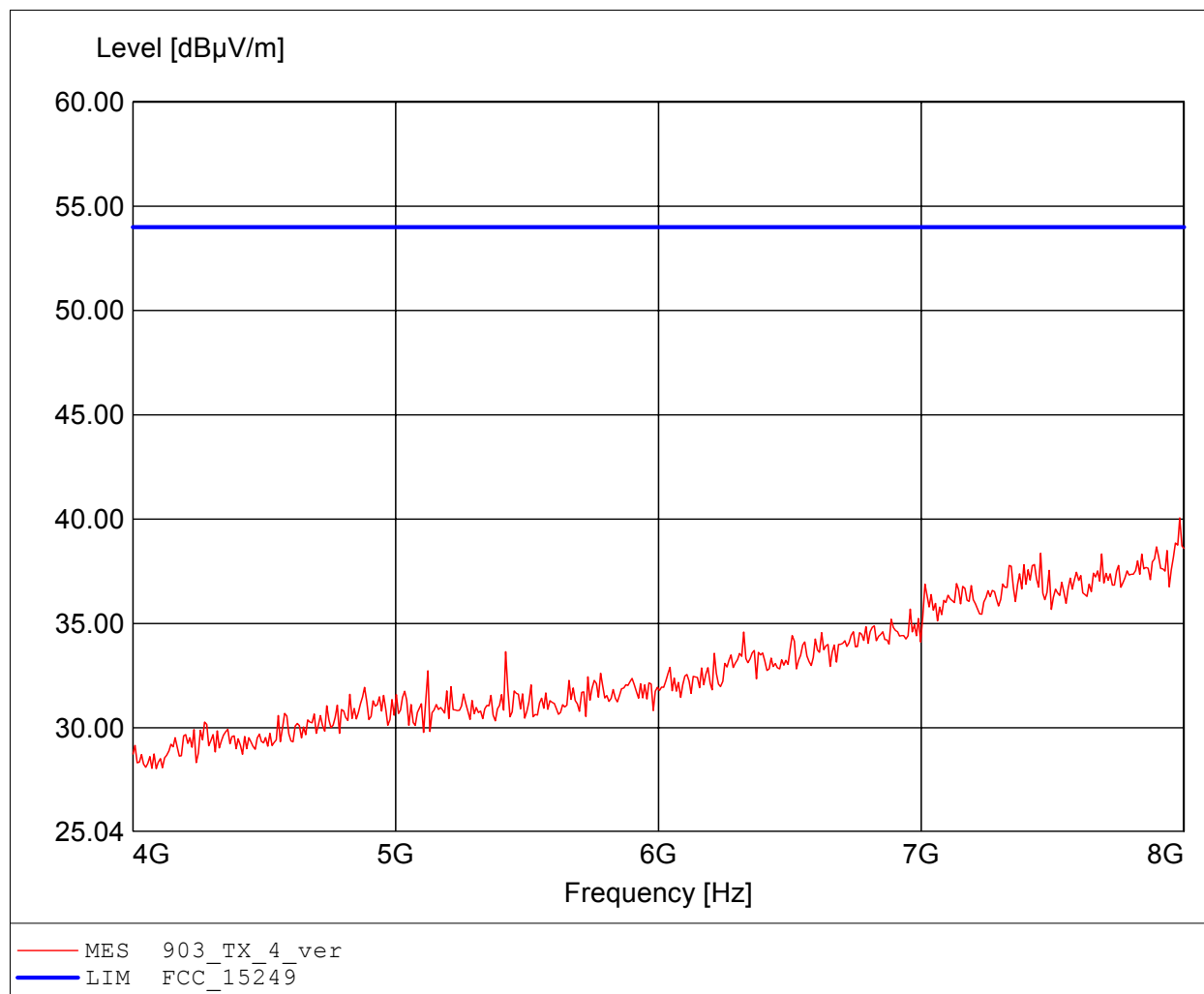

Spurious emissions Field Strength

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / G0M-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 1 / Freq: 903MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 903
Comment 1: Dist.: 3m, Ant.: BBHA9120D, ampl.+HP.
Comment 2: Freq: 5.419GHz, Emax: 41.18dBμV/m, RBW: 1MHz

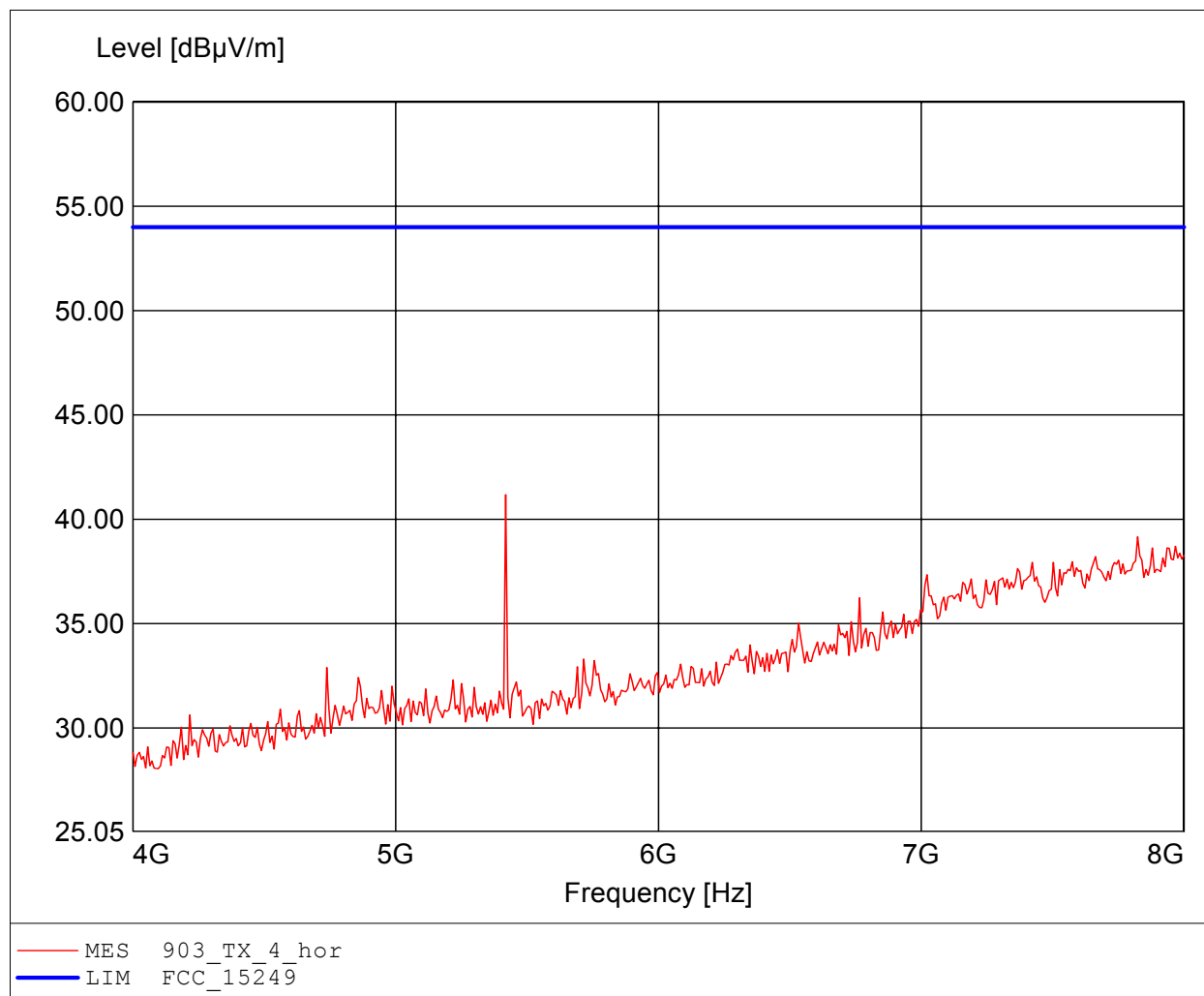

Spurious emissions Field Strength

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / GOM-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 1 / Freq: 903MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 903
Comment 1: Dist.: 3m, Ant.: BBHA9120D, ampl.+HP.
Comment 2: Freq: 11.944GHz, Emax: 42.57dBμV/m, RBW: 1MHz

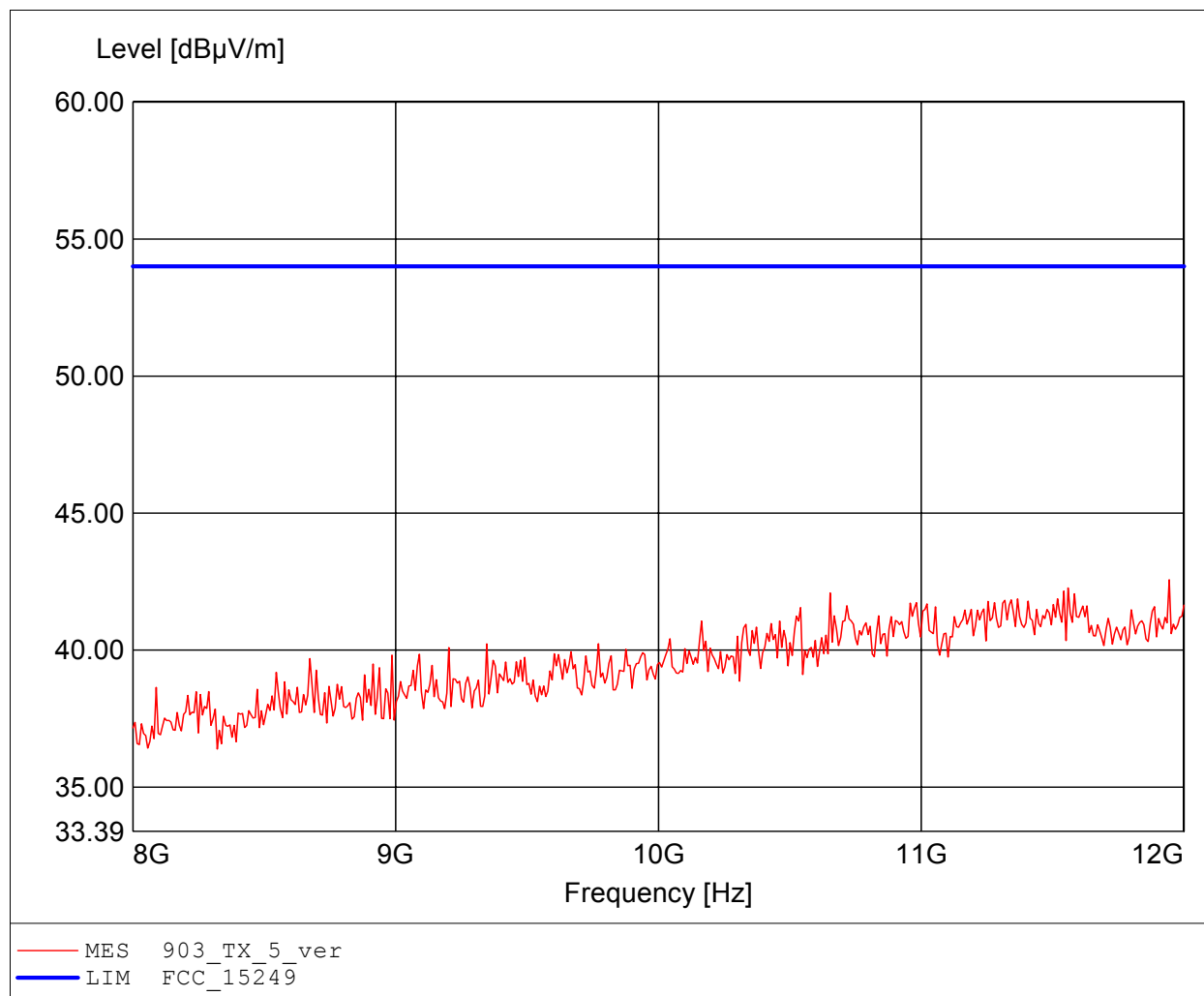

Spurious emissions Field Strength

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / G0M-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 1 / Freq: 903MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 903
Comment 1: Dist.: 3m, Ant.: BBHA9120D, ampl.+HP.
Comment 2: Freq: 11.287GHz, Emax: 42.80dBμV/m, RBW: 1MHz

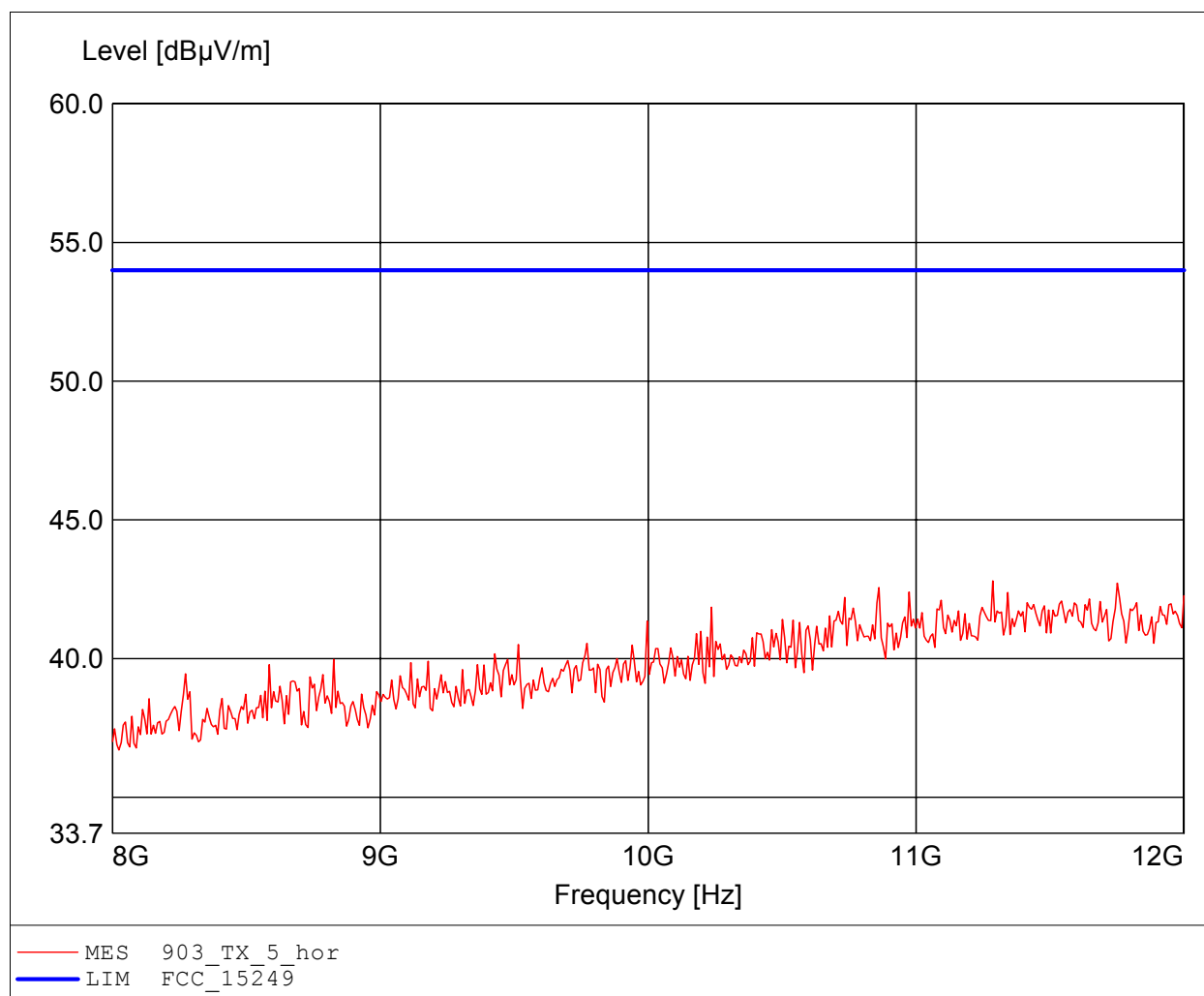

Spurious emissions Field Strength

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / GOM-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 8 / Freq: 911MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 911
Comment 1: Dist.: 3m, Ant.: HK 116
Comment 2: Freq: 188.758MHz, Emax: 33.06dBµV/m, RBW: 100kHz

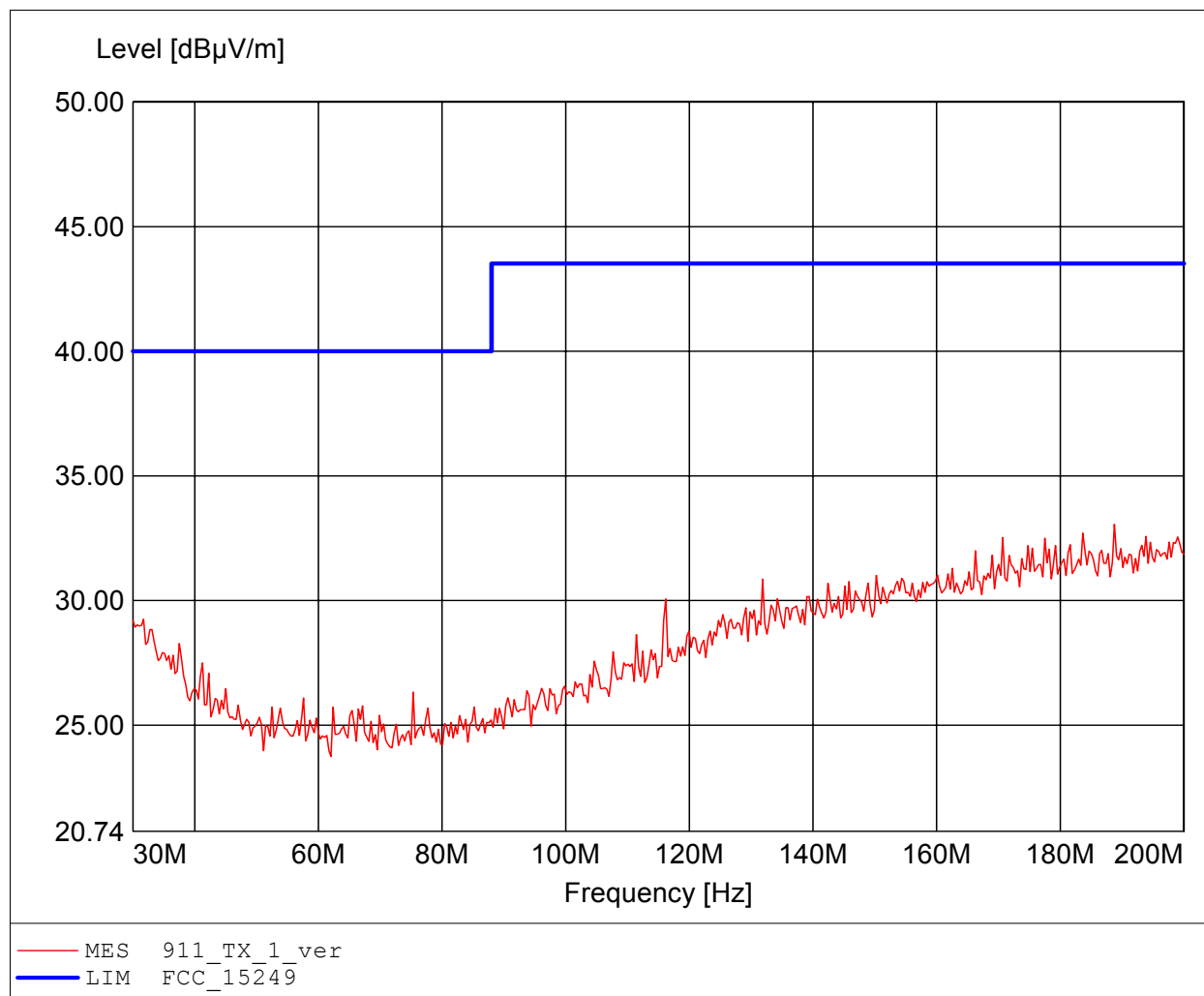

Spurious emissions Field Strength

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / G0M-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 8 / Freq: 911MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 911
Comment 1: Dist.: 3m, Ant.: HK 116
Comment 2: Freq: 189.098MHz, Emax: 33.19dBµV/m, RBW: 100kHz

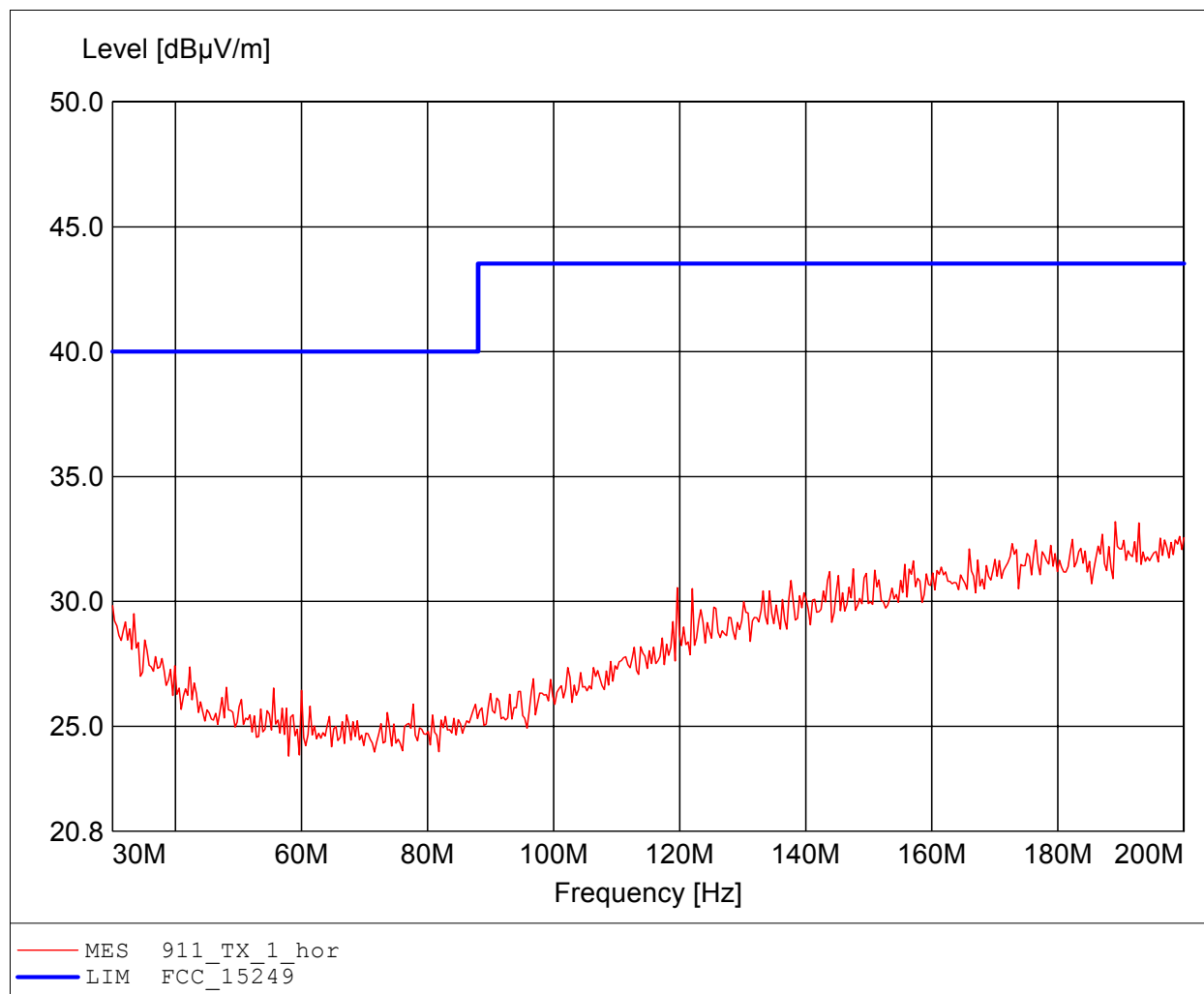

Spurious emissions Field Strength

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / GOM-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 8 / Freq: 911MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 911
Comment 1: Dist.: 3m, Ant.: HL 223, amplif.
Comment 2: Freq: 902.000MHz, Emax: 26.30dBµV/m, RBW: 100kHz

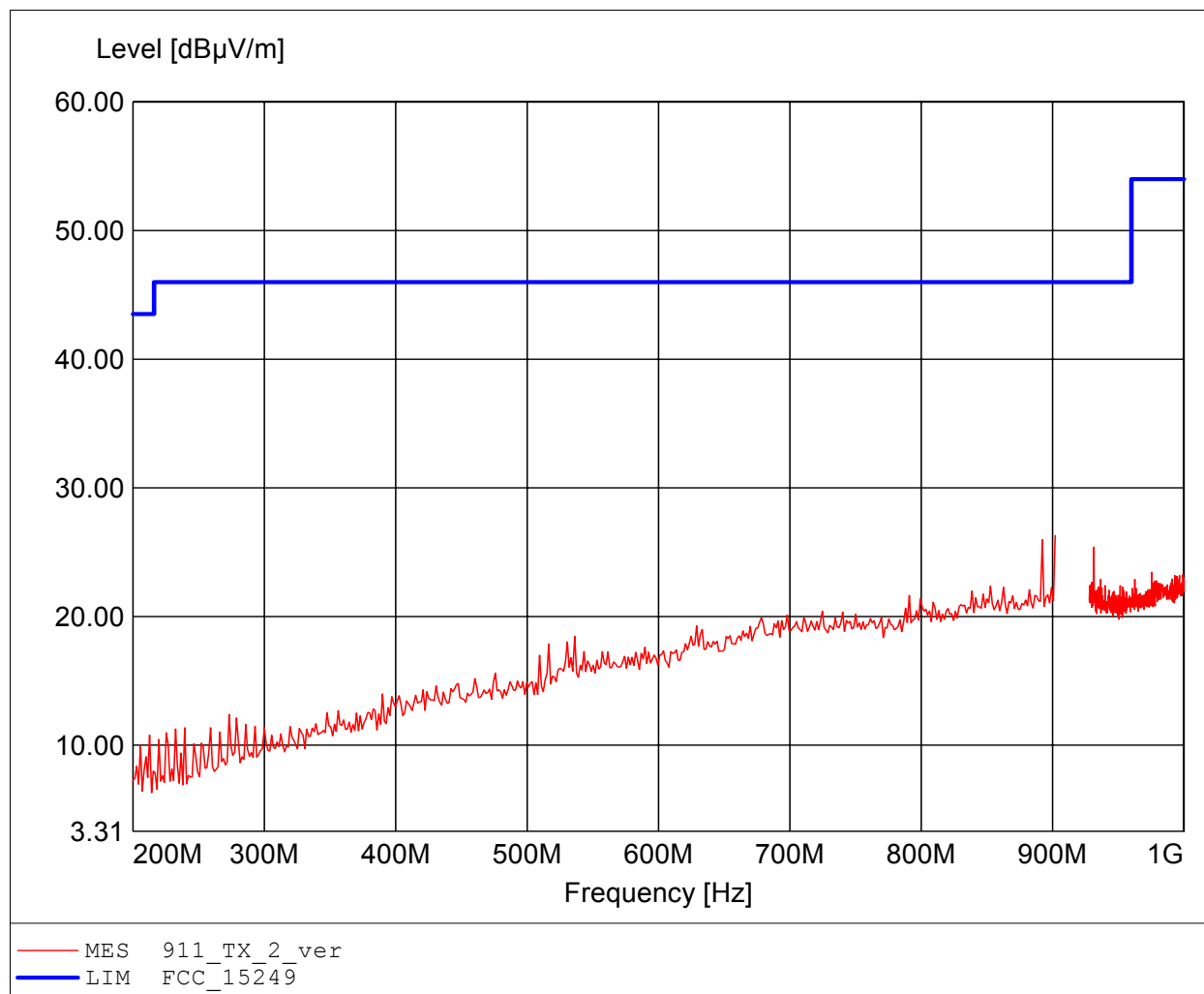

Spurious emissions Field Strength

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / G0M-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 8 / Freq: 911MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 911
Comment 1: Dist.: 3m, Ant.: HL 223, amplif.
Comment 2: Freq: 931.463MHz, Emax: 28.15dBµV/m, RBW: 100kHz

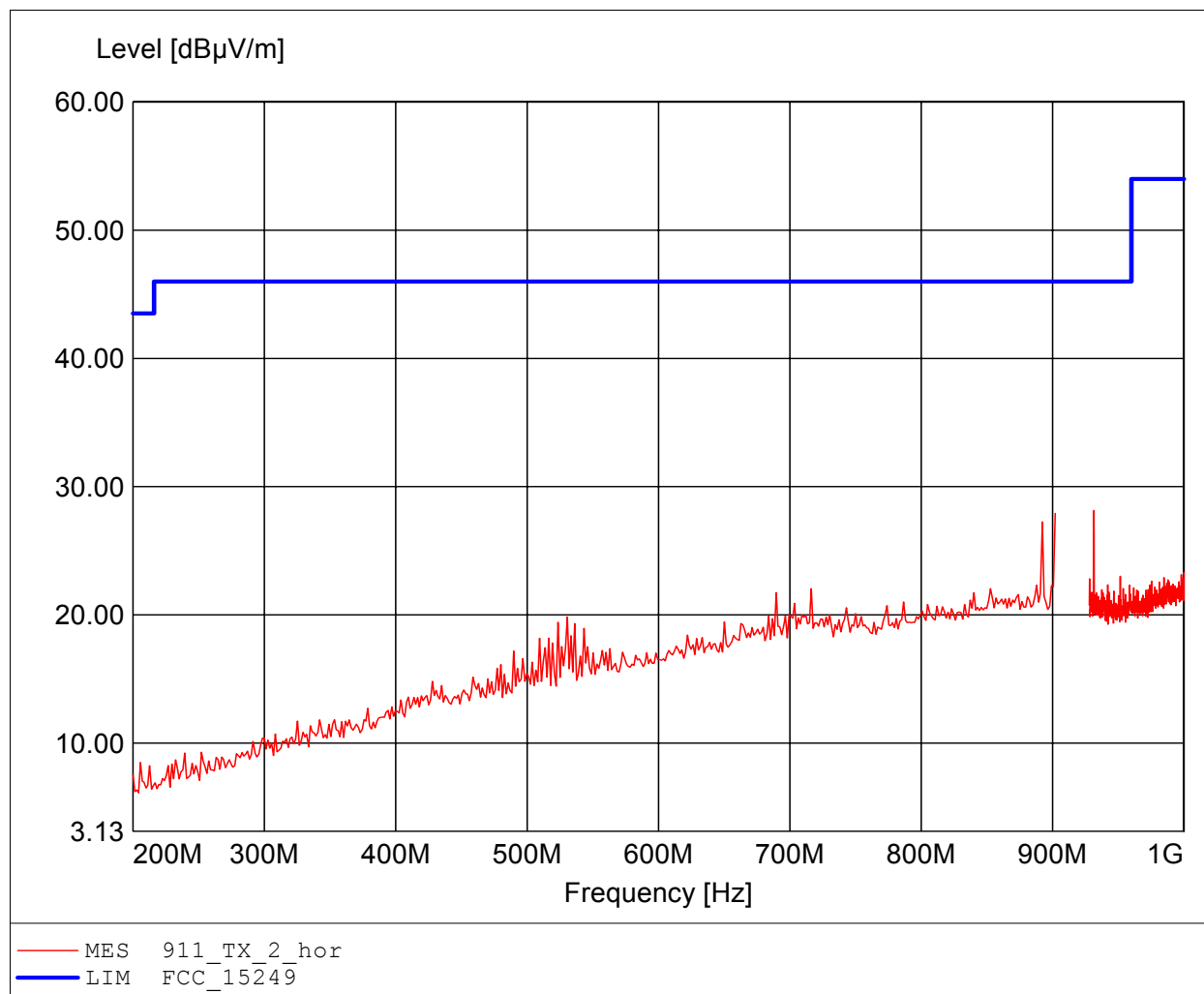

Spurious emissions Field Strength

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / GOM-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 8 / Freq: 911MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 911
Comment 1: Dist.: 3m, Ant.: BBHA9120D, amplif.
Comment 2: Freq: 3.916GHz, Emax: 42.96dB μ V/m, RBW: 1MHz

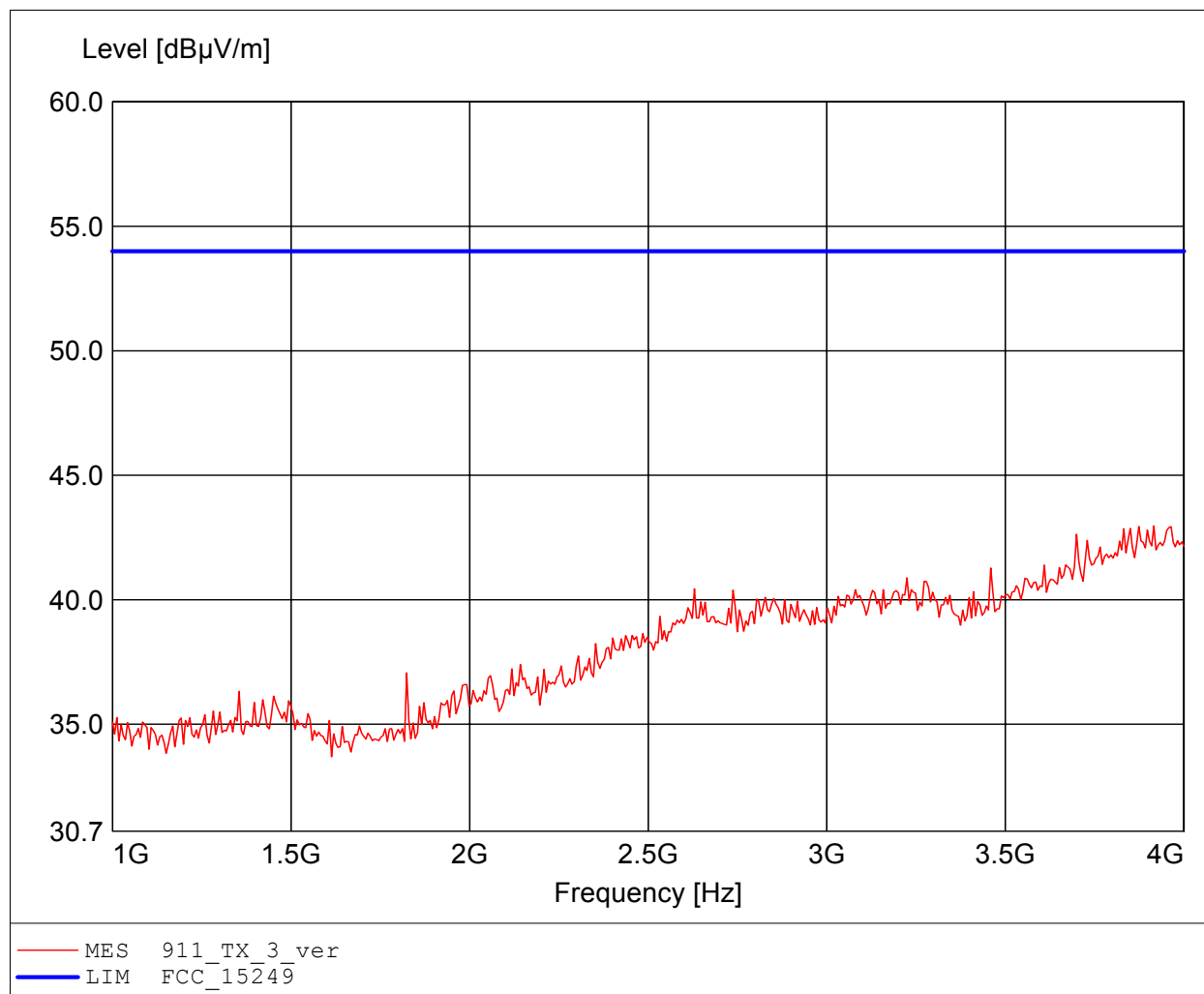

Spurious emissions Field Strength

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / G0M-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 8 / Freq: 911MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 911
Comment 1: Dist.: 3m, Ant.: BBHA9120D, amplif.
Comment 2: Freq: 3.850GHz, Emax: 43.77dBμV/m, RBW: 1MHz

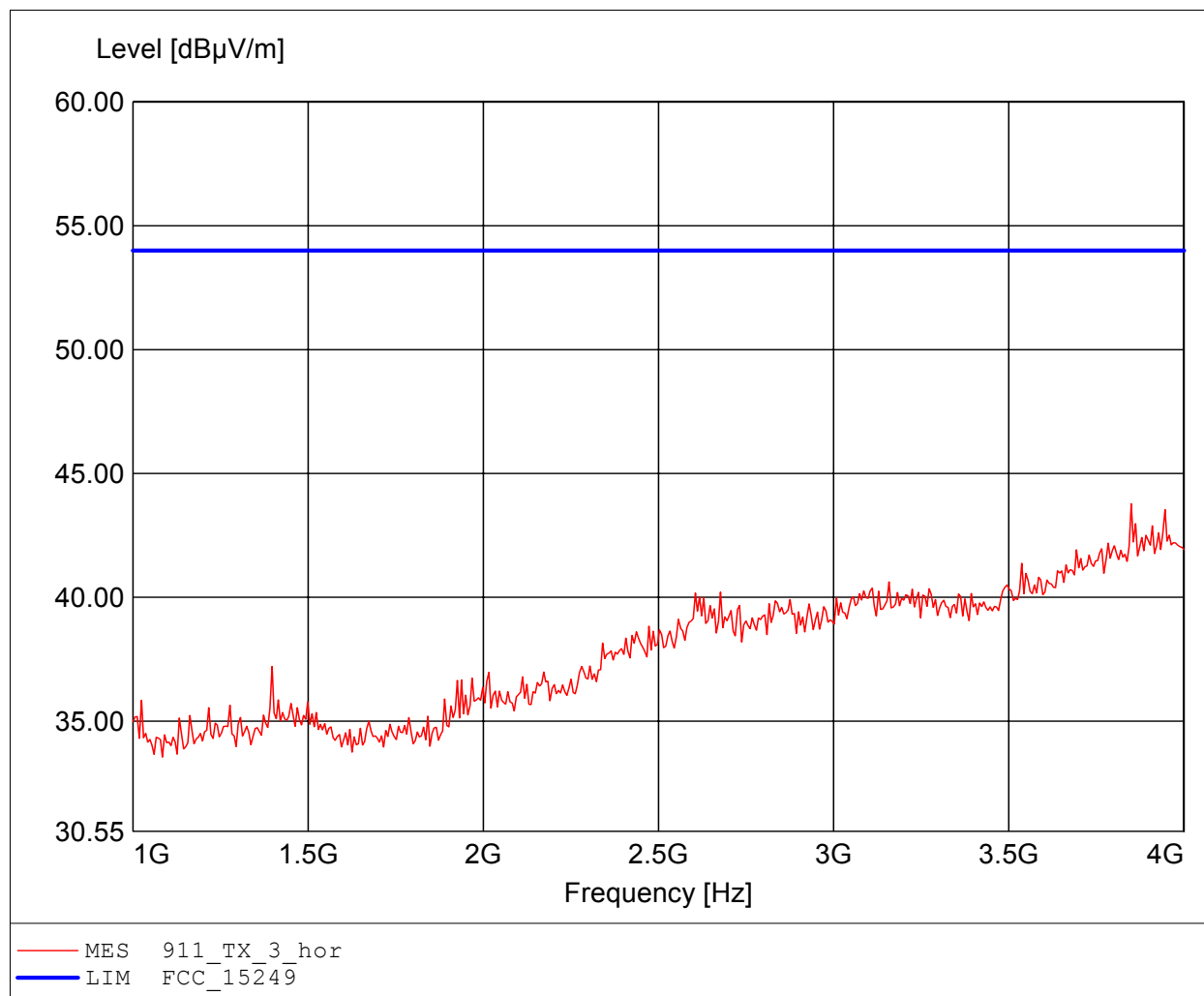

Spurious emissions Field Strength

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / GOM-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 8 / Freq: 911MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 911
Comment 1: Dist.: 3m, Ant.: BBHA9120D, ampl.+HP.
Comment 2: Freq: 7.992GHz, Emax: 39.31dBµV/m, RBW: 1MHz

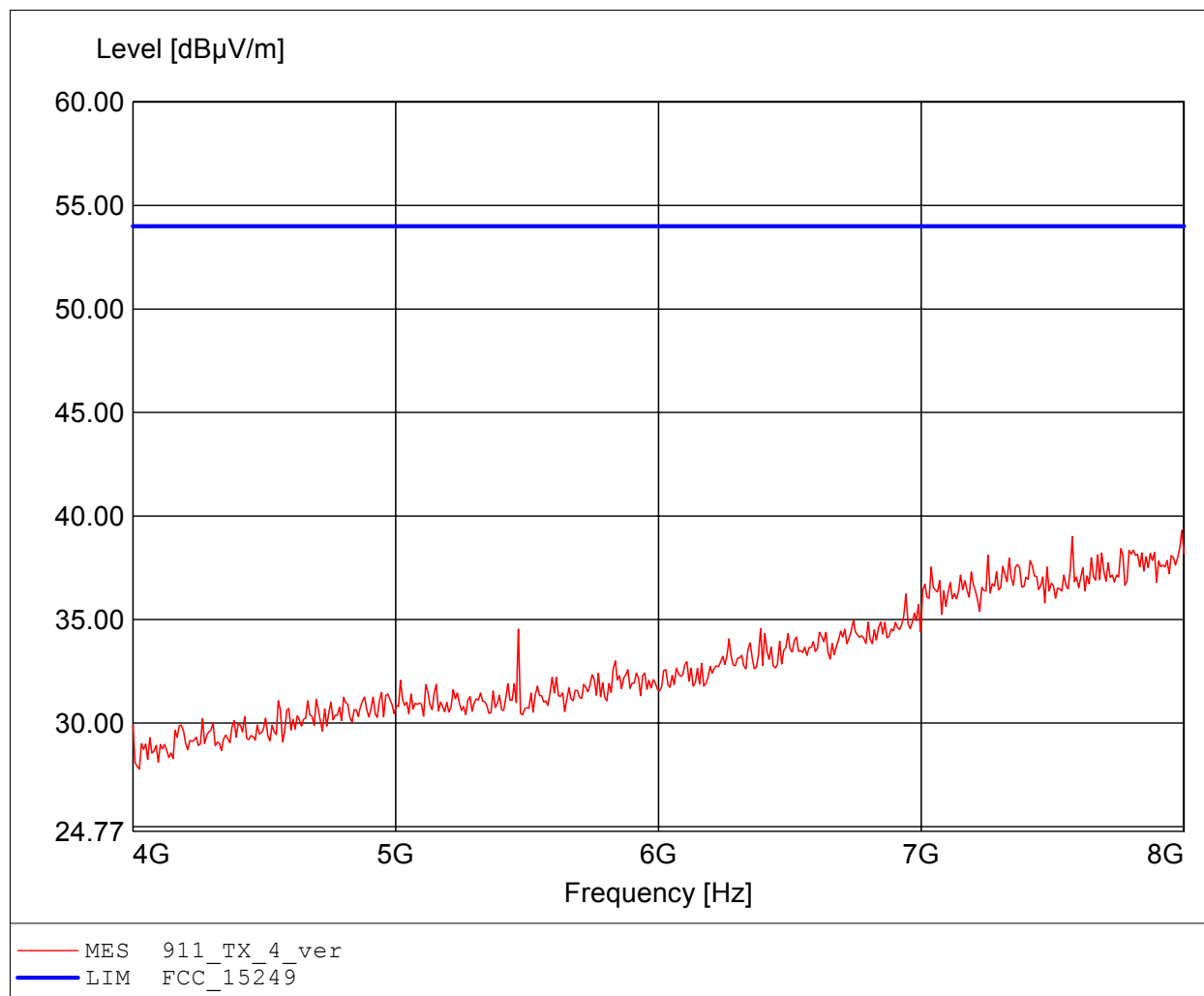

Spurious emissions Field Strength

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / G0M-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 8 / Freq: 911MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 911
Comment 1: Dist.: 3m, Ant.: BBHA9120D, ampl.+HP.
Comment 2: Freq: 5.467GHz, Emax: 41.44dBµV/m, RBW: 1MHz

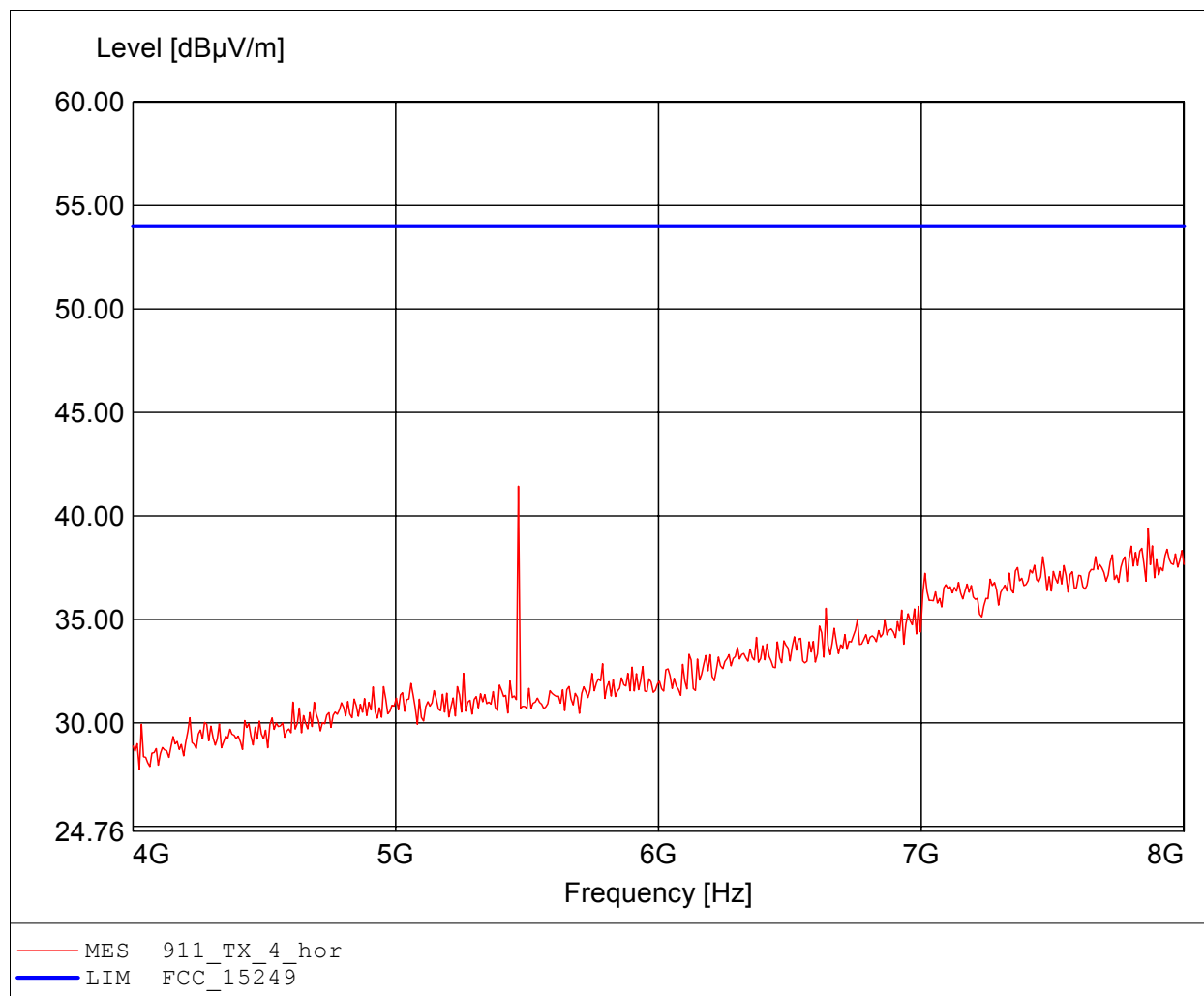

Spurious emissions Field Strength

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / GOM-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 8 / Freq: 911MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 911
Comment 1: Dist.: 3m, Ant.: BBHA9120D, ampl.+HP.
Comment 2: Freq: 11.776GHz, Emax: 43.24dB μ V/m, RBW: 1MHz

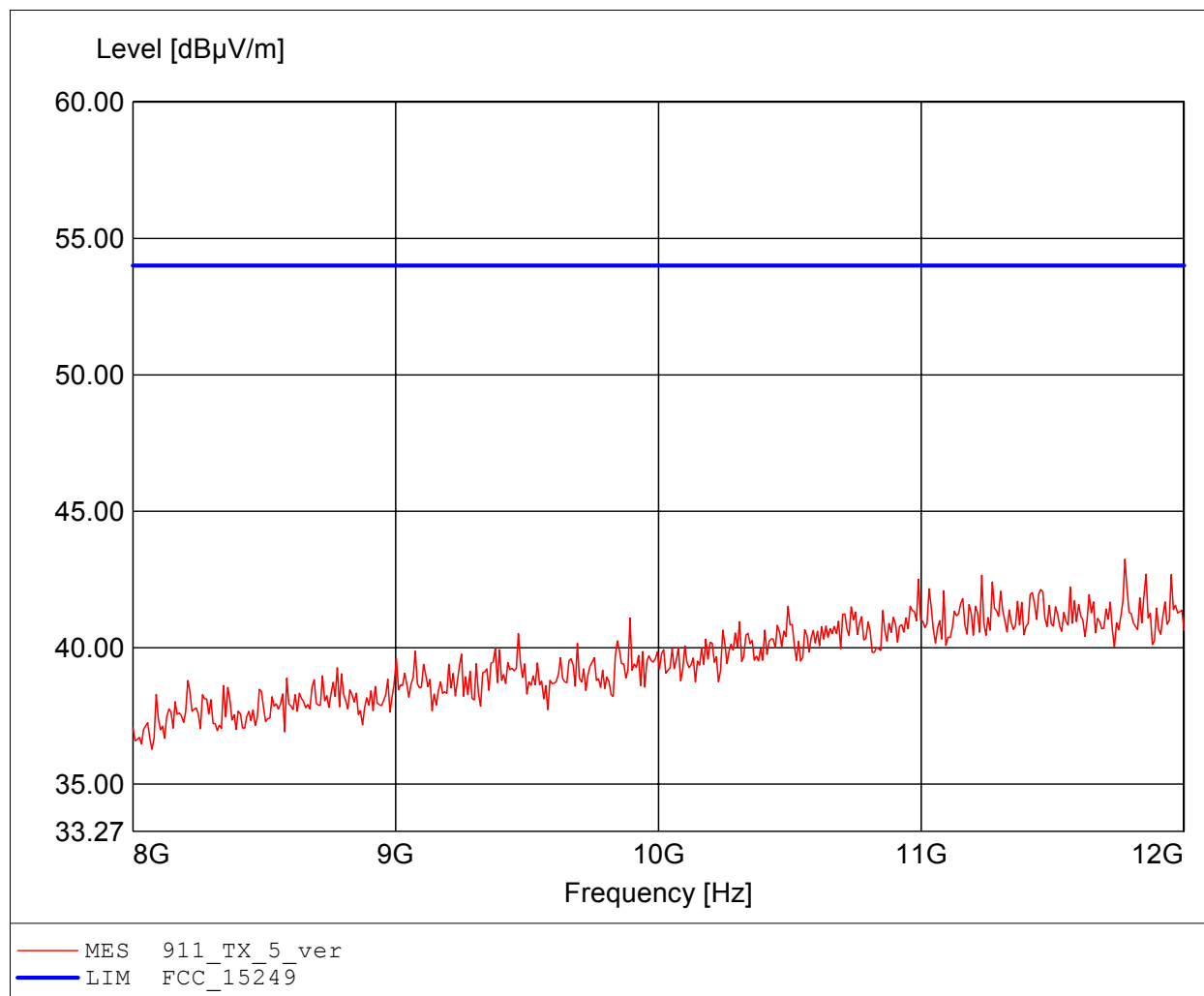

Spurious emissions Field Strength

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / G0M-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 8 / Freq: 911MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 911
Comment 1: Dist.: 3m, Ant.: BBHA9120D, ampl.+HP.
Comment 2: Freq: 11.359GHz, Emax: 42.67dBµV/m, RBW: 1MHz

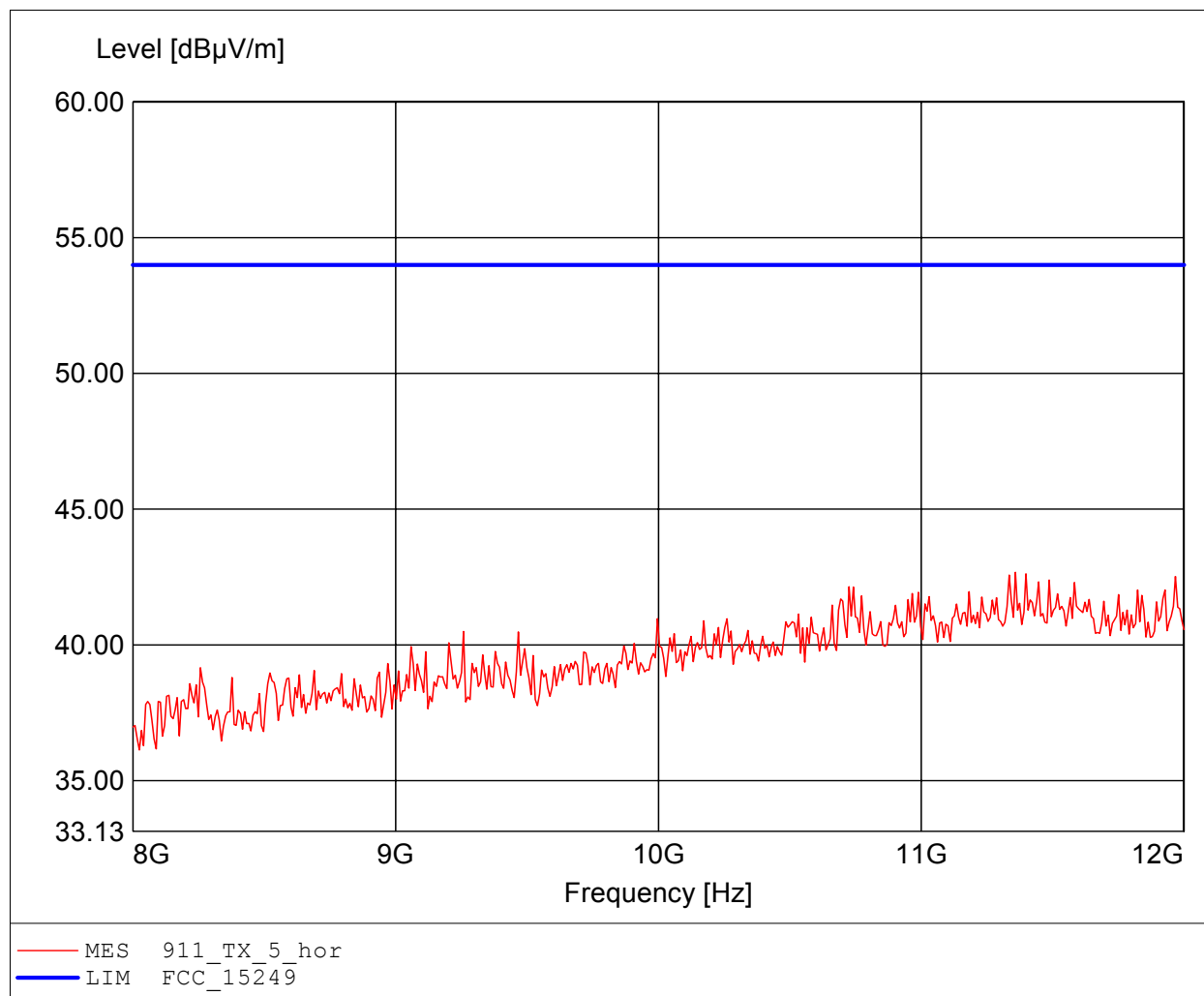

Spurious emissions Field Strength

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / GOM-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 16 / Freq: 921MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 921
Comment 1: Dist.: 3m, Ant.: HK 116
Comment 2: Freq: 196.253MHz, Emax: 32.75dBµV/m, RBW: 100kHz

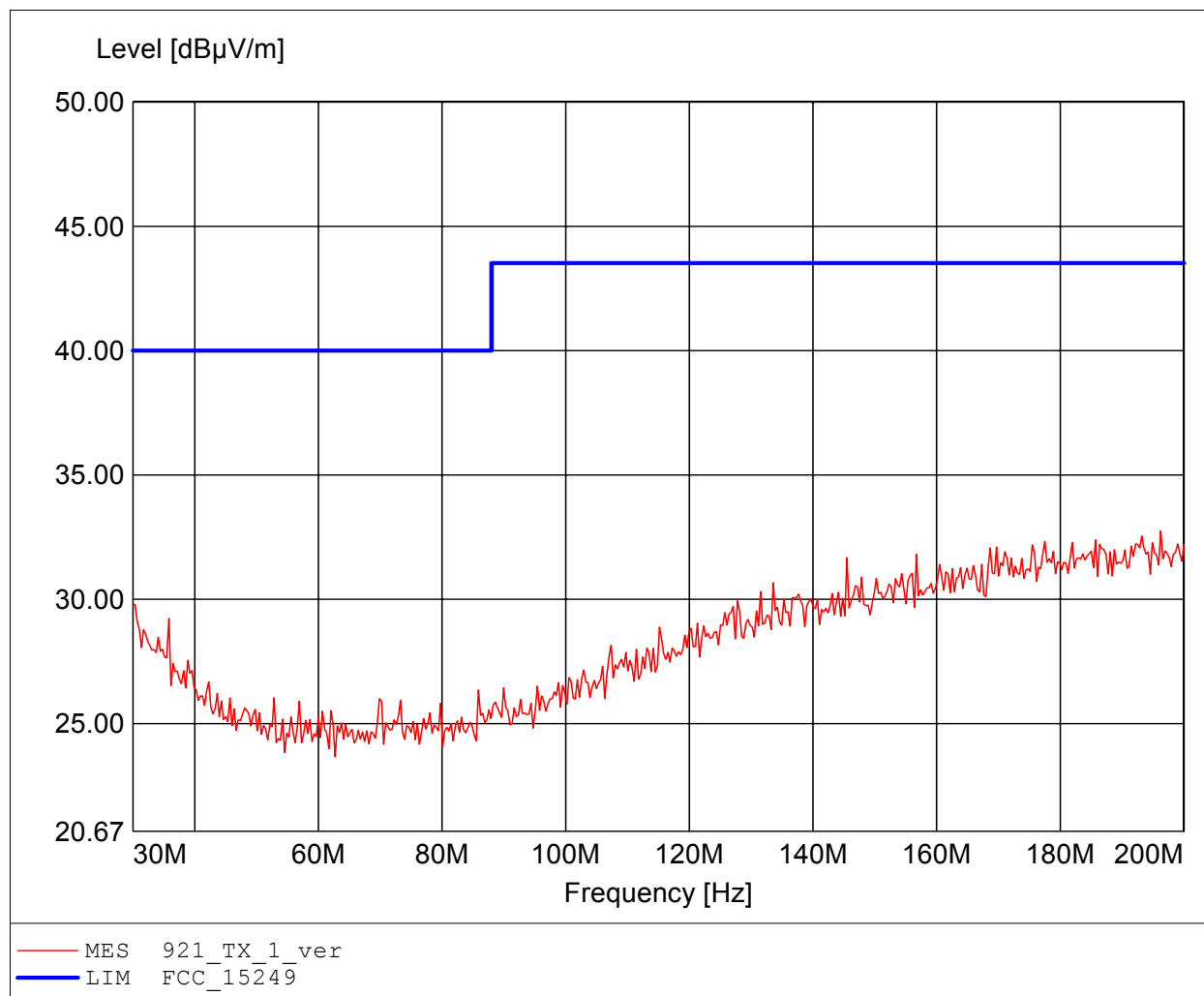

Spurious emissions Field Strength

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / G0M-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 16 / Freq: 921MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 921
Comment 1: Dist.: 3m, Ant.: HK 116
Comment 2: Freq: 193.868MHz, Emax: 33.06dBµV/m, RBW: 100kHz

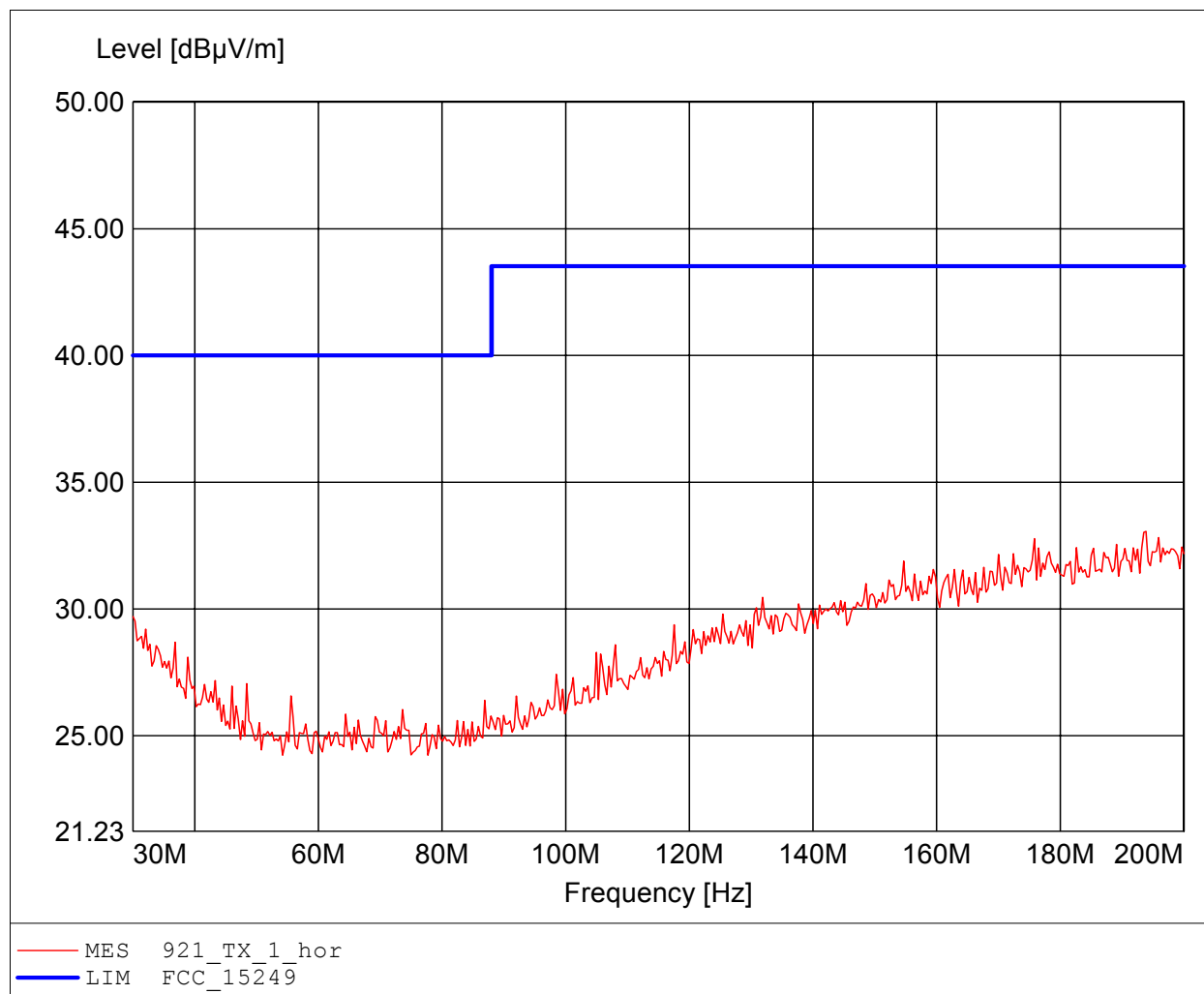

Spurious emissions Field Strength

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / GOM-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 16 / Freq: 921MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 921
Comment 1: Dist.: 3m, Ant.: HL 223, amplif.
Comment 2: Freq: 931.174MHz, Emax: 30.89dBµV/m, RBW: 100kHz

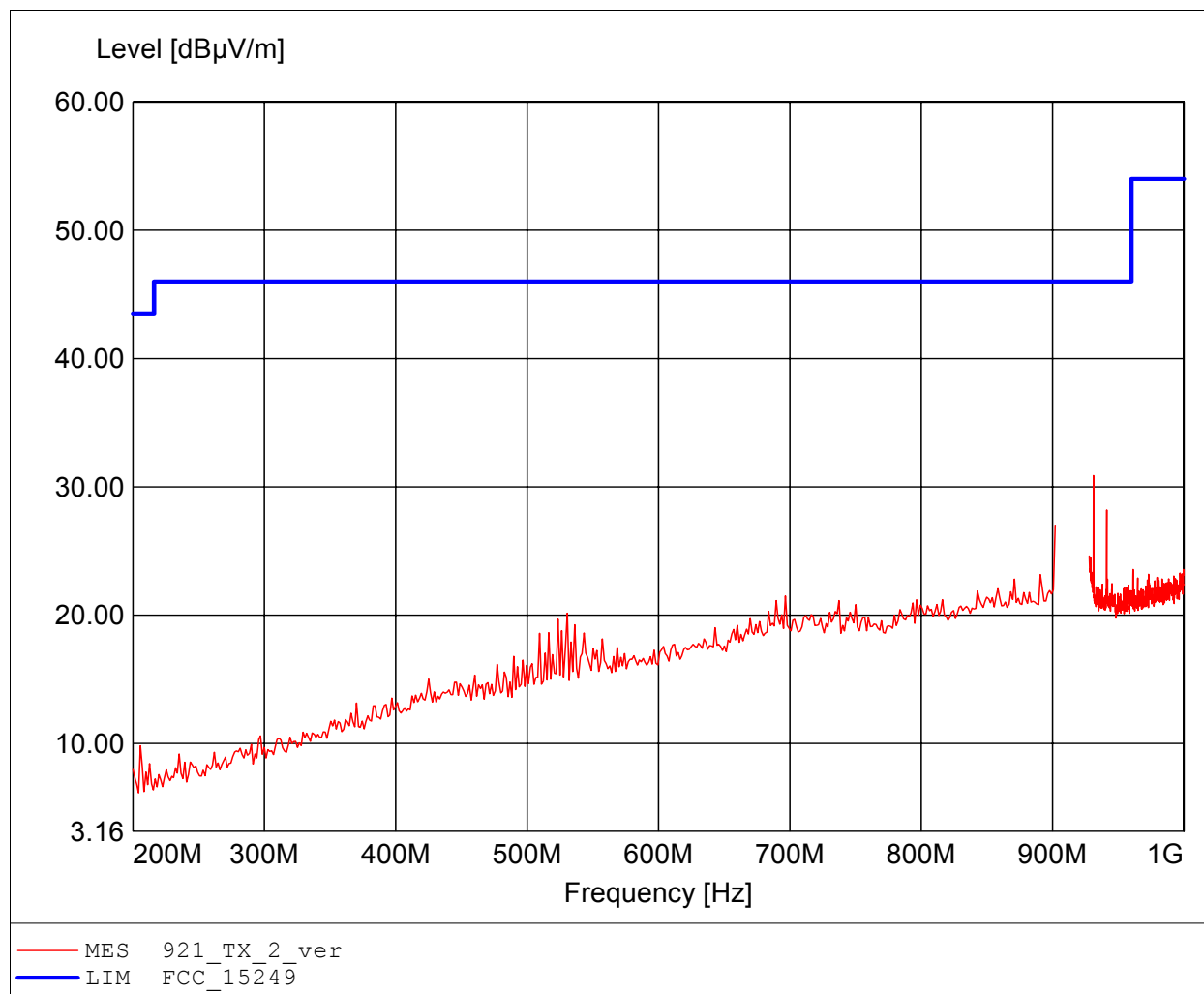

Spurious emissions Field Strength

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / G0M-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 16 / Freq: 921MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 921
Comment 1: Dist.: 3m, Ant.: HL 223, amplif.
Comment 2: Freq: 931.174MHz, Emax: 26.71dBµV/m, RBW: 100kHz

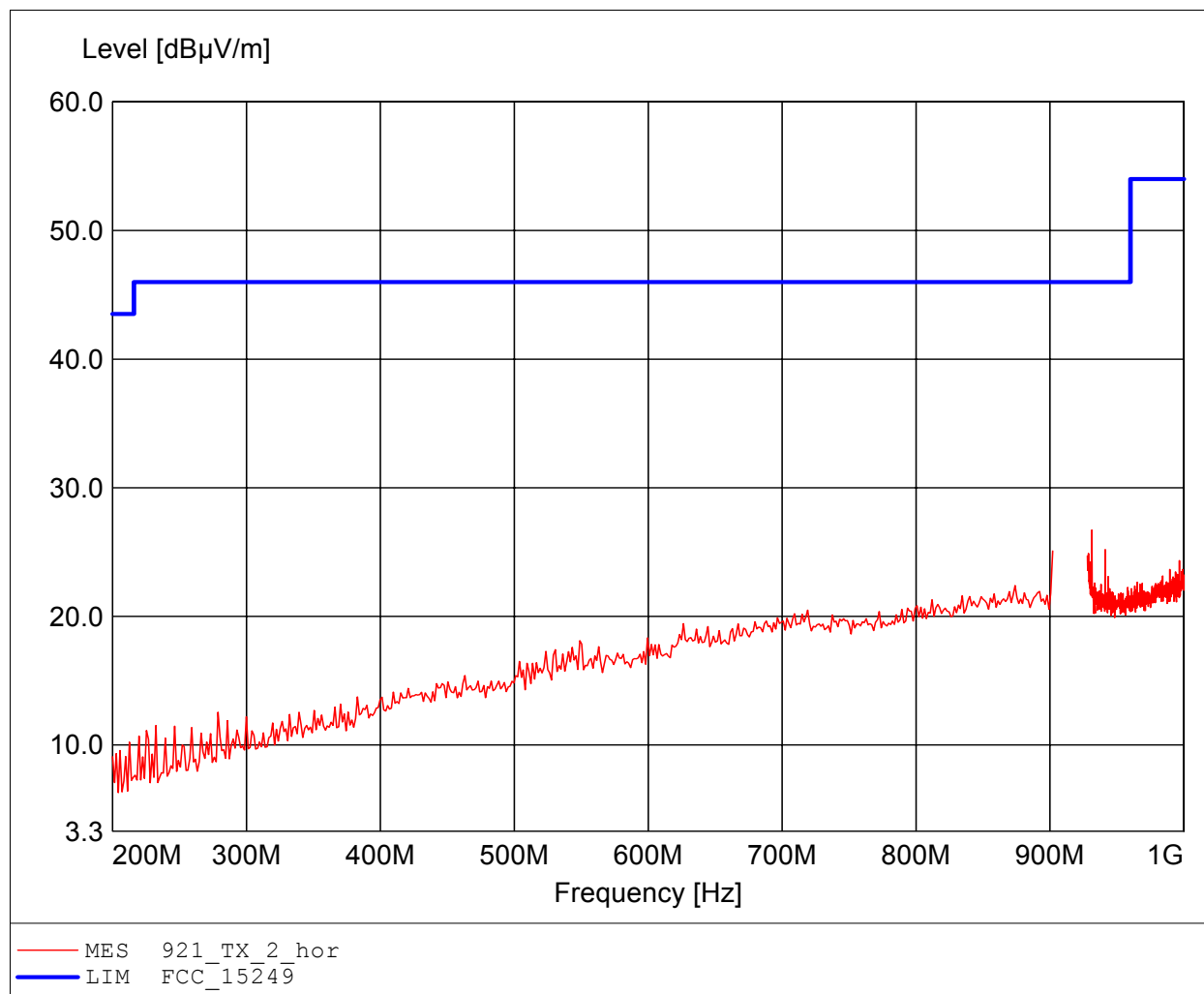

Spurious emissions Field Strength

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / GOM-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 16 / Freq: 921MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 921
Comment 1: Dist.: 3m, Ant.: BBHA9120D, amplif.
Comment 2: Freq: 3.826GHz, Emax: 42.57dBμV/m, RBW: 1MHz

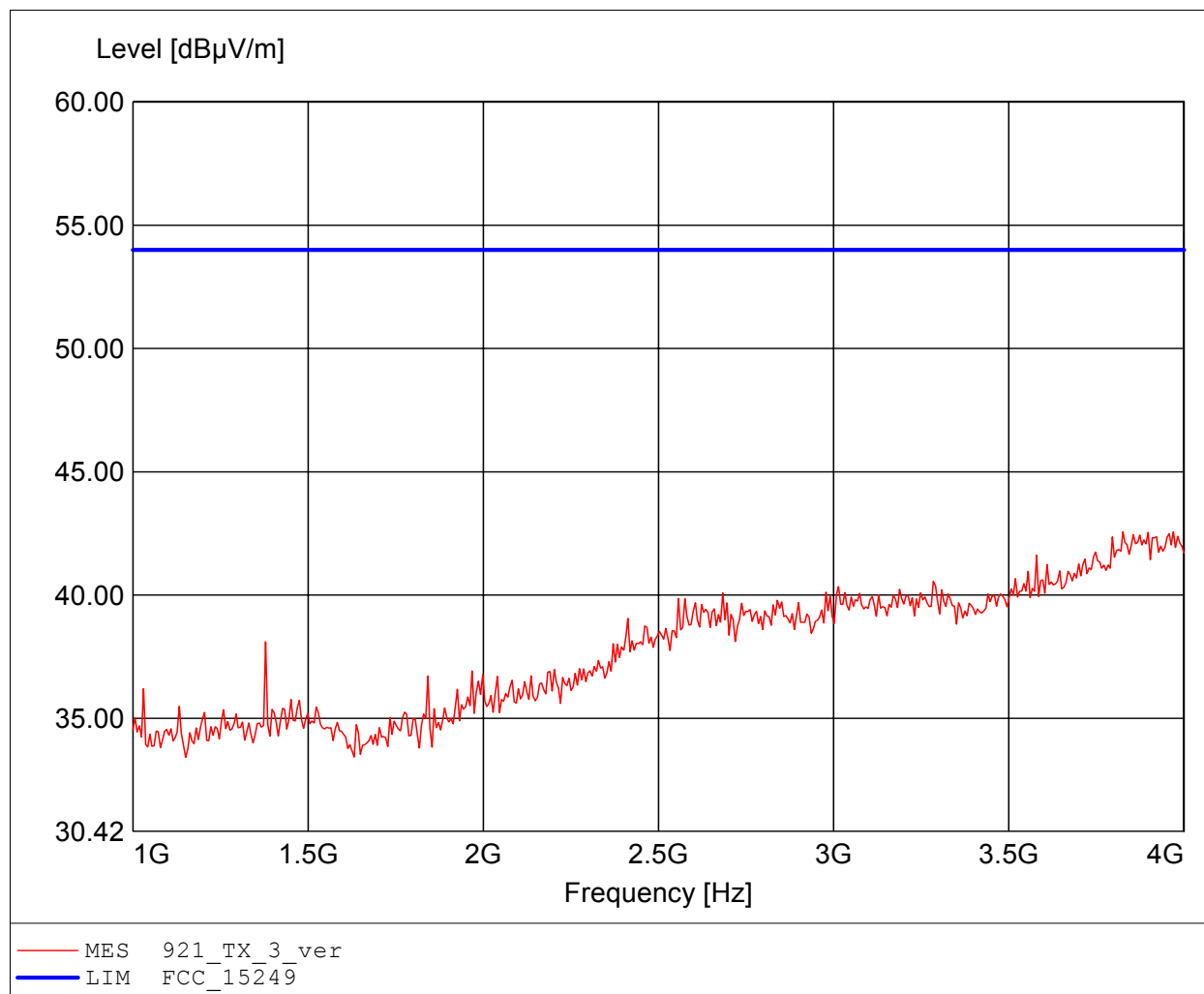

Spurious emissions Field Strength

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / G0M-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 16 / Freq: 921MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 921
Comment 1: Dist.: 3m, Ant.: BBHA9120D, amplif.
Comment 2: Freq: 3.946GHz, Emax: 42.62dBµV/m, RBW: 1MHz

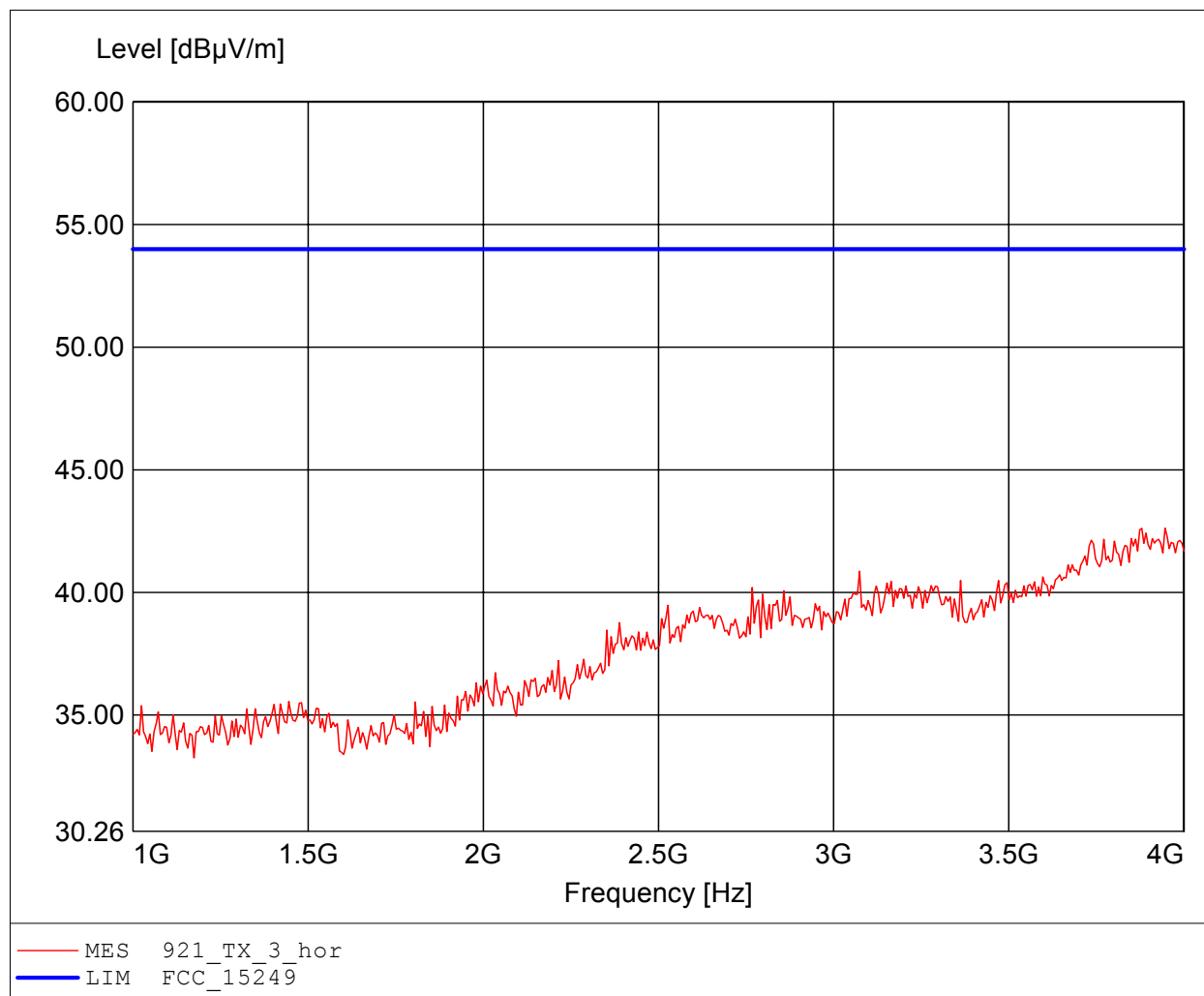

Spurious emissions Field Strength

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / GOM-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 16 / Freq: 921MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 921
Comment 1: Dist.: 3m, Ant.: BBHA9120D, ampl.+HP.
Comment 2: Freq: 7.984GHz, Emax: 39.32dBµV/m, RBW: 1MHz

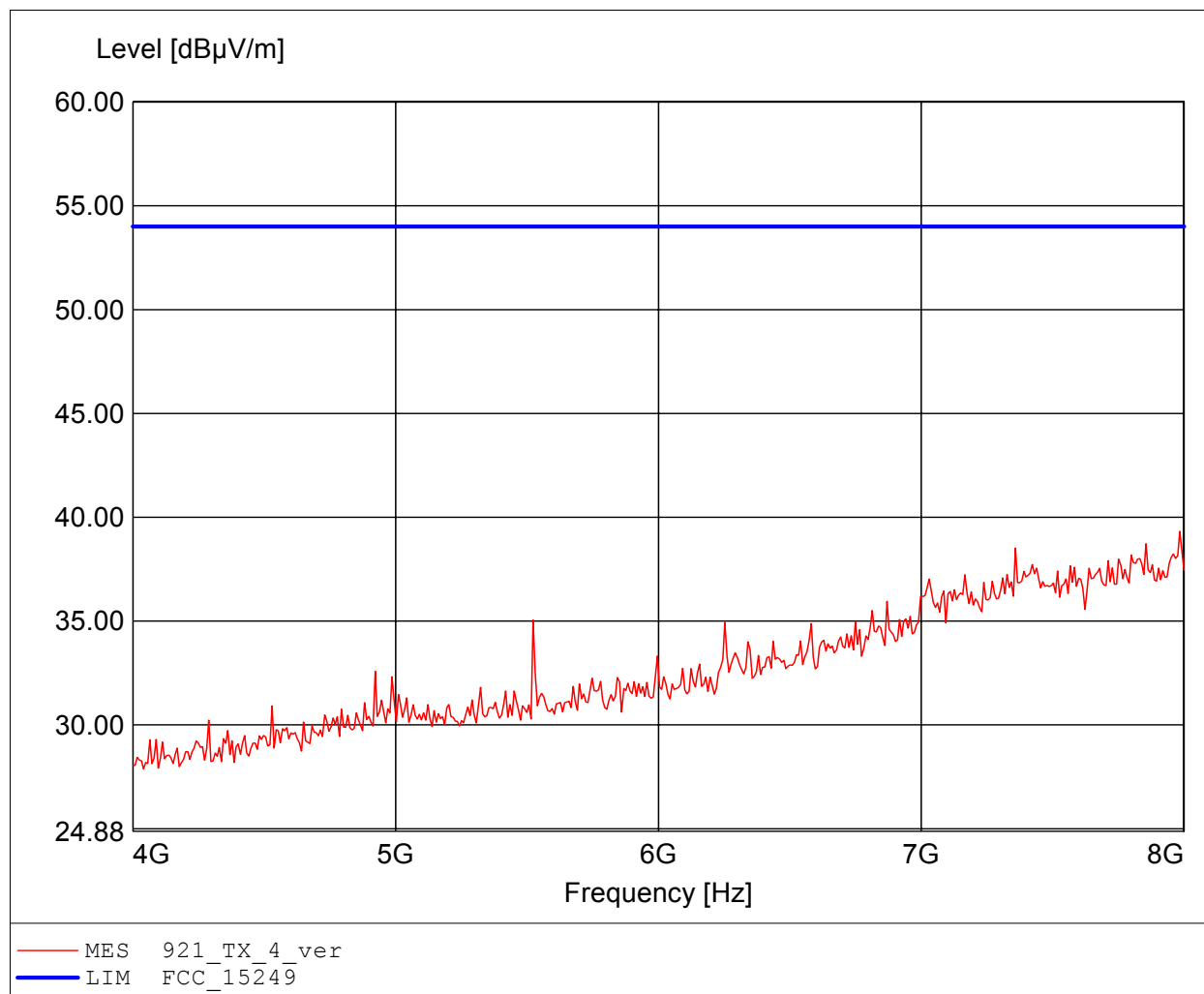

Spurious emissions Field Strength

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / G0M-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 16 / Freq: 921MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 921
Comment 1: Dist.: 3m, Ant.: BBHA9120D, ampl.+HP.
Comment 2: Freq: 5.523GHz, Emax: 41.63dBμV/m, RBW: 1MHz

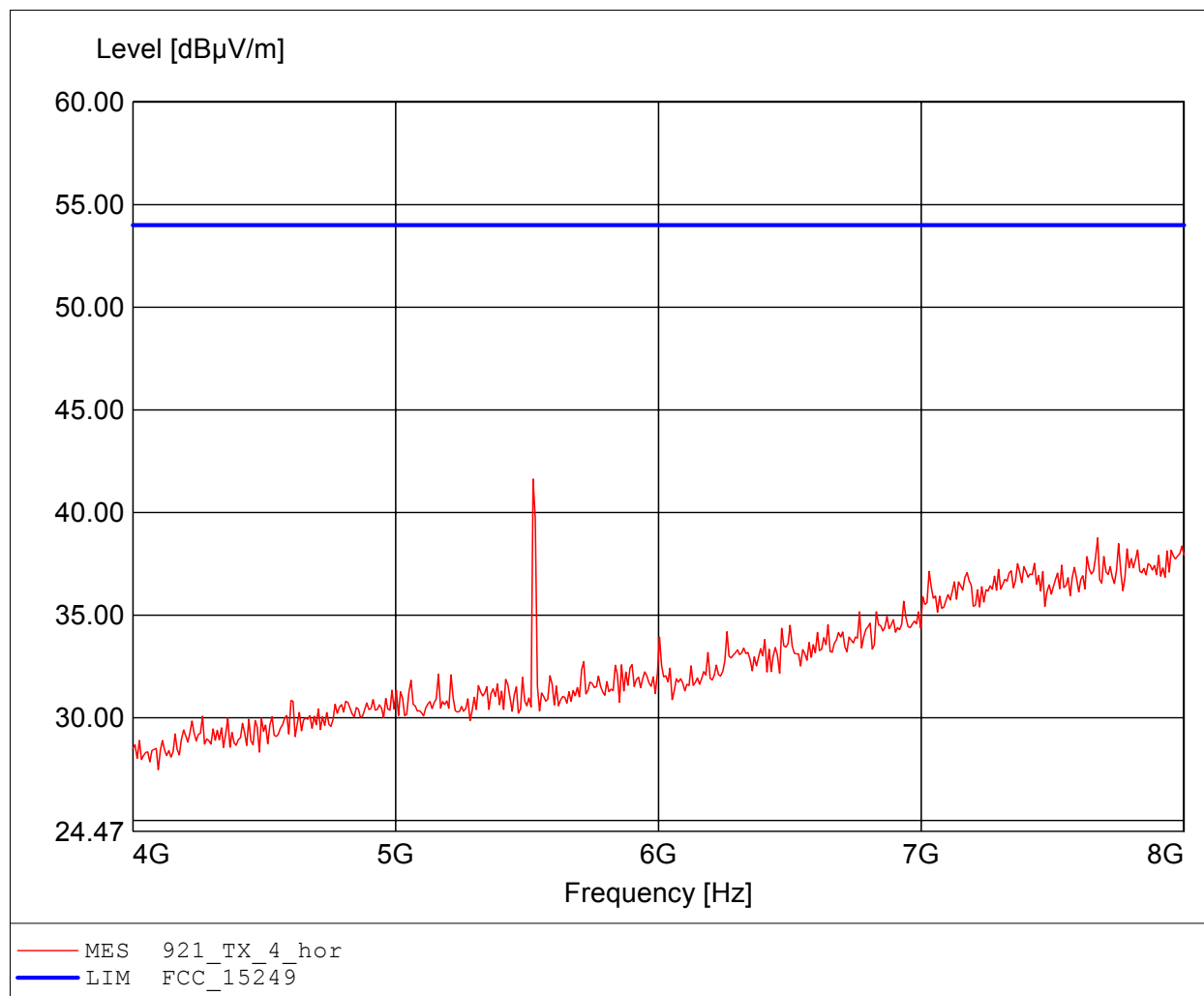

Spurious emissions Field Strength

FCC RULES PART 15, SUBPART C


Approval Holder: LifeScan Scotland Ltd. / GOM-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 16 / Freq: 921MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 921
Comment 1: Dist.: 3m, Ant.: BBHA9120D, ampl.+HP.
Comment 2: Freq: 11.543GHz, Emax: 43.05dBμV/m, RBW: 1MHz

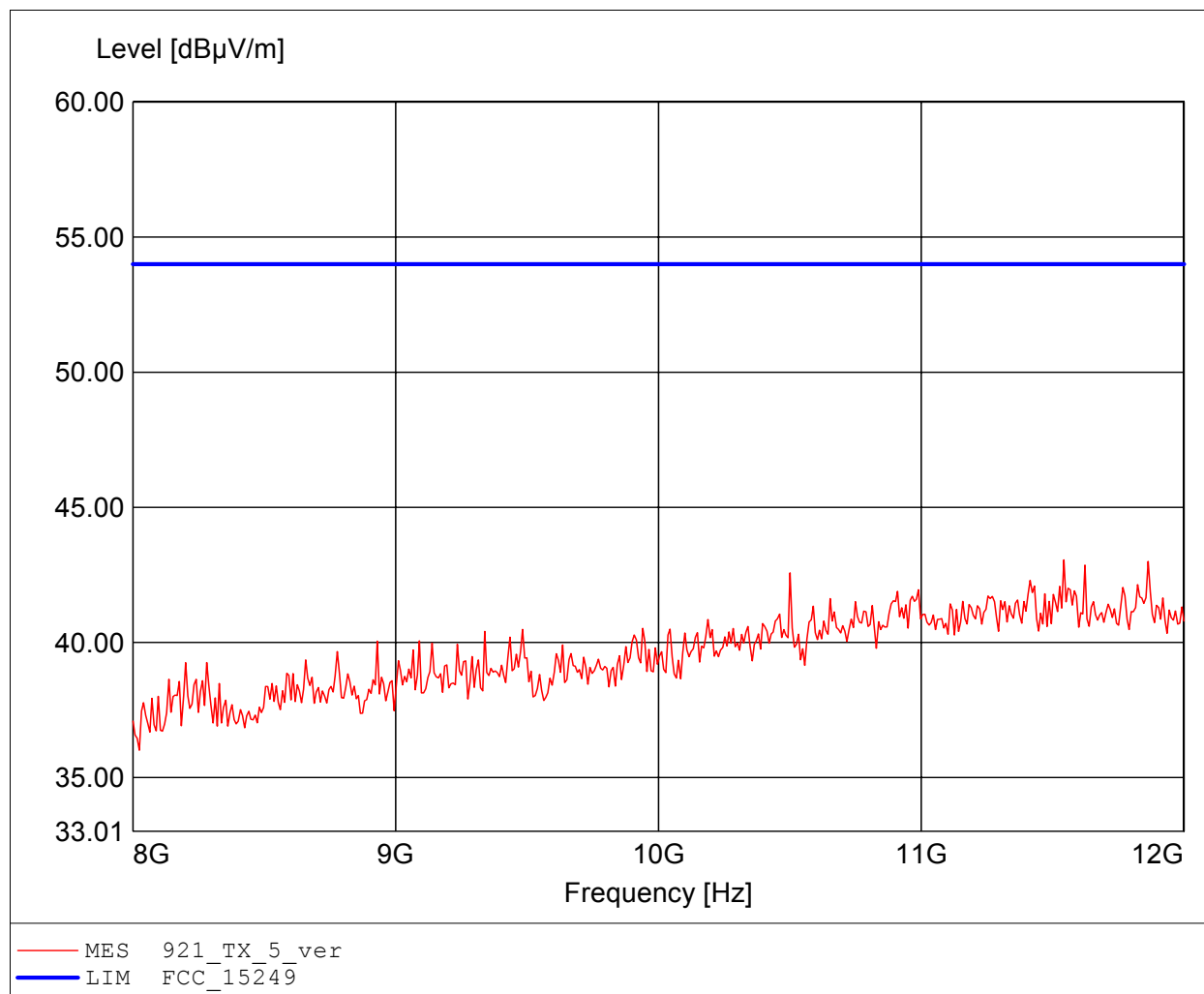

Spurious emissions Field Strength

FCC RULES PART 15, SUBPART C

Approval Holder: LifeScan Scotland Ltd. / G0M-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 16 / Freq: 921MHz / FSK / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 921
Comment 1: Dist.: 3m, Ant.: BBHA9120D, ampl.+HP.
Comment 2: Freq: 11.447GHz, Emax: 43.05dBμV/m, RBW: 1MHz

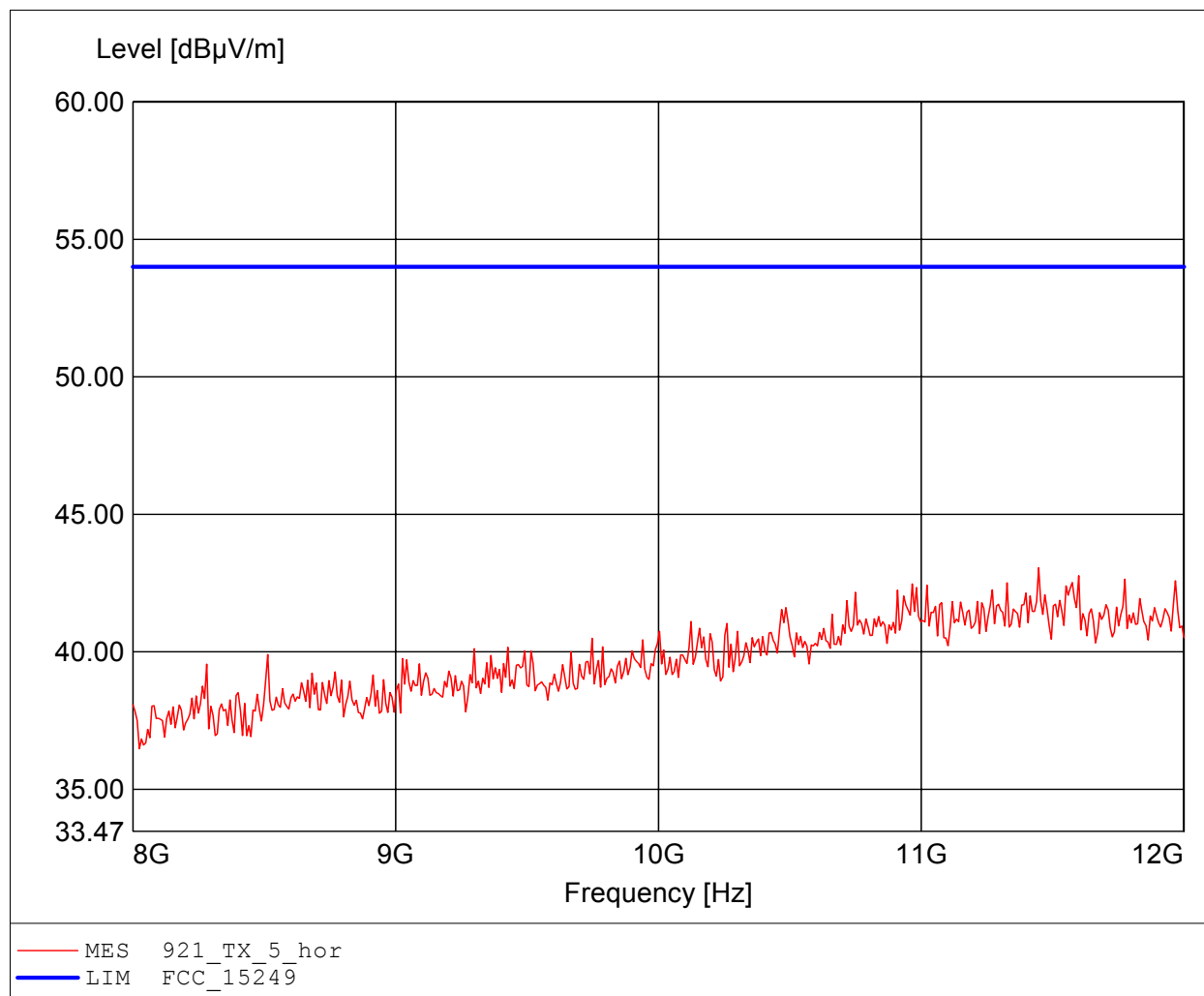


Annex E Receiver radiated spurious emissions

Field Strength under normal conditions

Standards Industry Canada, RSS-GEN


Approval Holder: LifeScan Scotland Ltd. / G0M-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 8 / Freq: 911MHz / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 8
Comment 1: Dist.: 3m, Ant.: HK 116
Comment 2: Freq:199.659MHz Emax:33.14dBµV/m RBW: 100 kHz

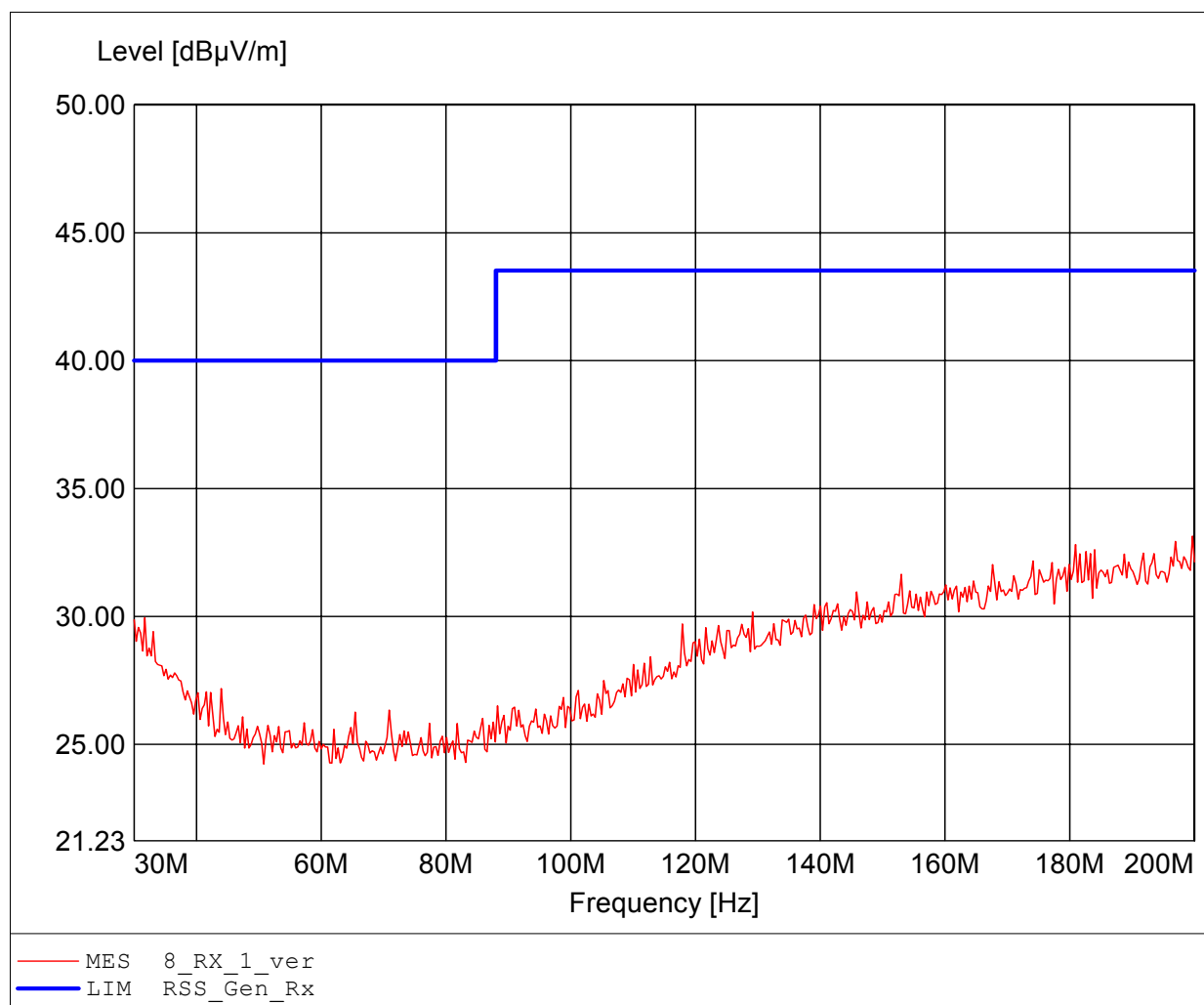

Field Strength under normal conditions

Standards Industry Canada, RSS-GEN


Approval Holder: LifeScan Scotland Ltd. / G0M-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 8 / Freq: 911MHz / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 8
Comment 1: Dist.: 3m, Ant.: HK 116
Comment 2: Freq:191.142MHz Emax:33.19dBµV/m RBW: 100 kHz

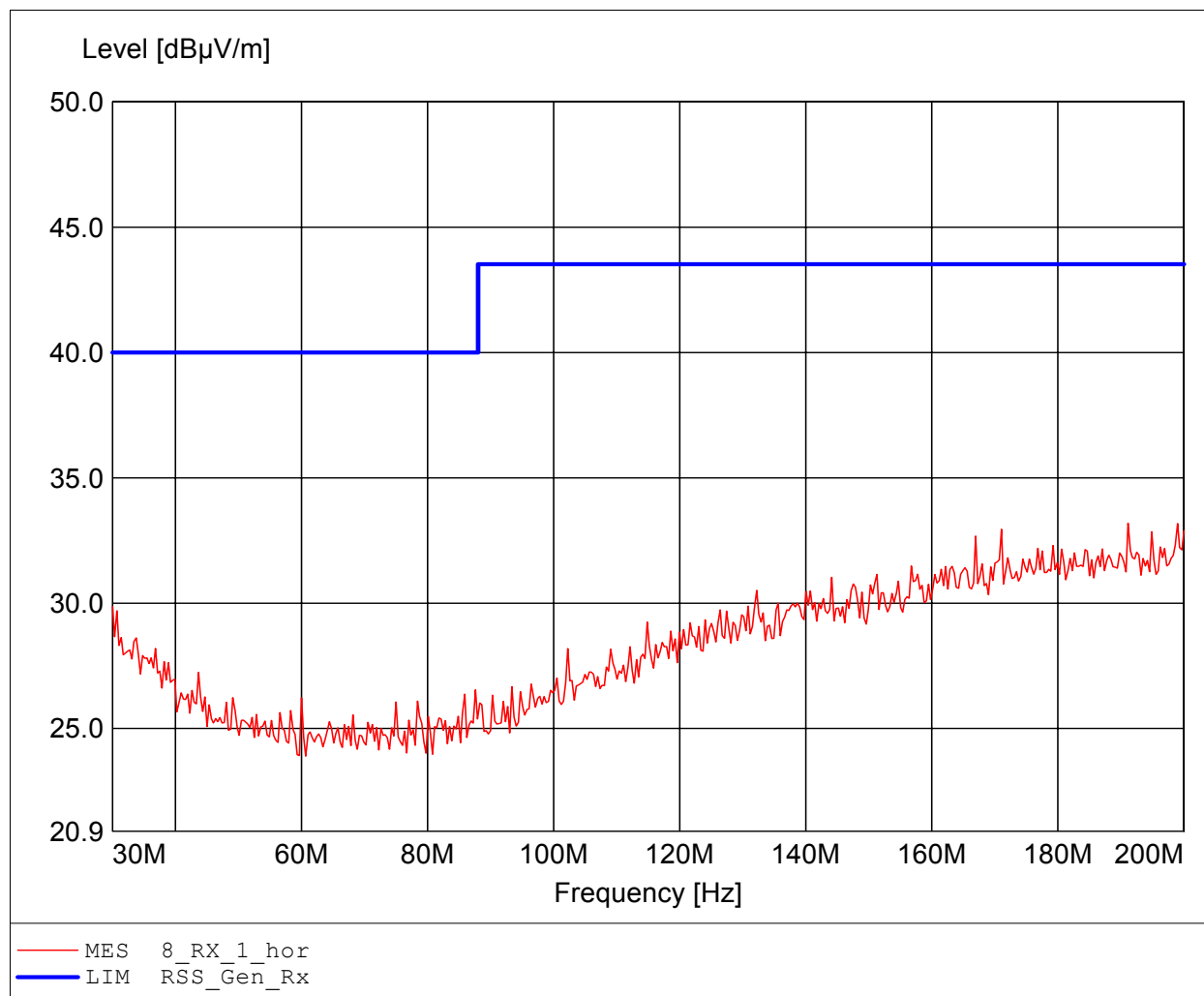

Field Strength under normal conditions

Standards Industry Canada, RSS-GEN


Approval Holder: LifeScan Scotland Ltd. / G0M-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 8 / Freq: 911MHz / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 8
Comment 1: Dist.: 3m, Ant.: HL 223, ampl.
Comment 2: Freq: 897.395MHz Emax: 34.04dBµV/m RBW: 100 kHz

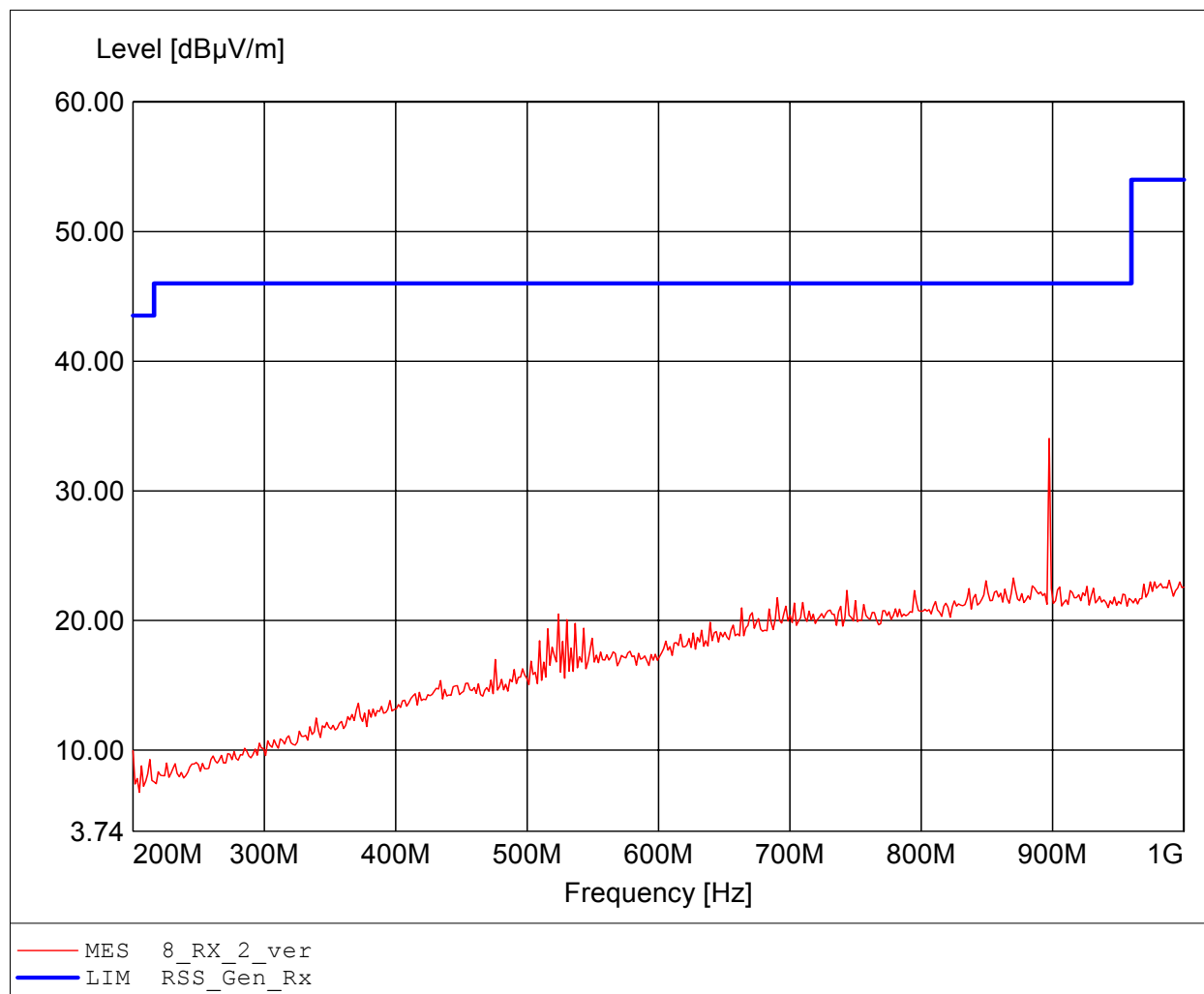

Field Strength under normal conditions

Standards Industry Canada, RSS-GEN


Approval Holder: LifeScan Scotland Ltd. / G0M-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 8 / Freq: 911MHz / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 8
Comment 1: Dist.: 3m, Ant.: HL 223, ampl.
Comment 2: Freq: 898.998MHz Emax: 30.24dBµV/m RBW: 100 kHz

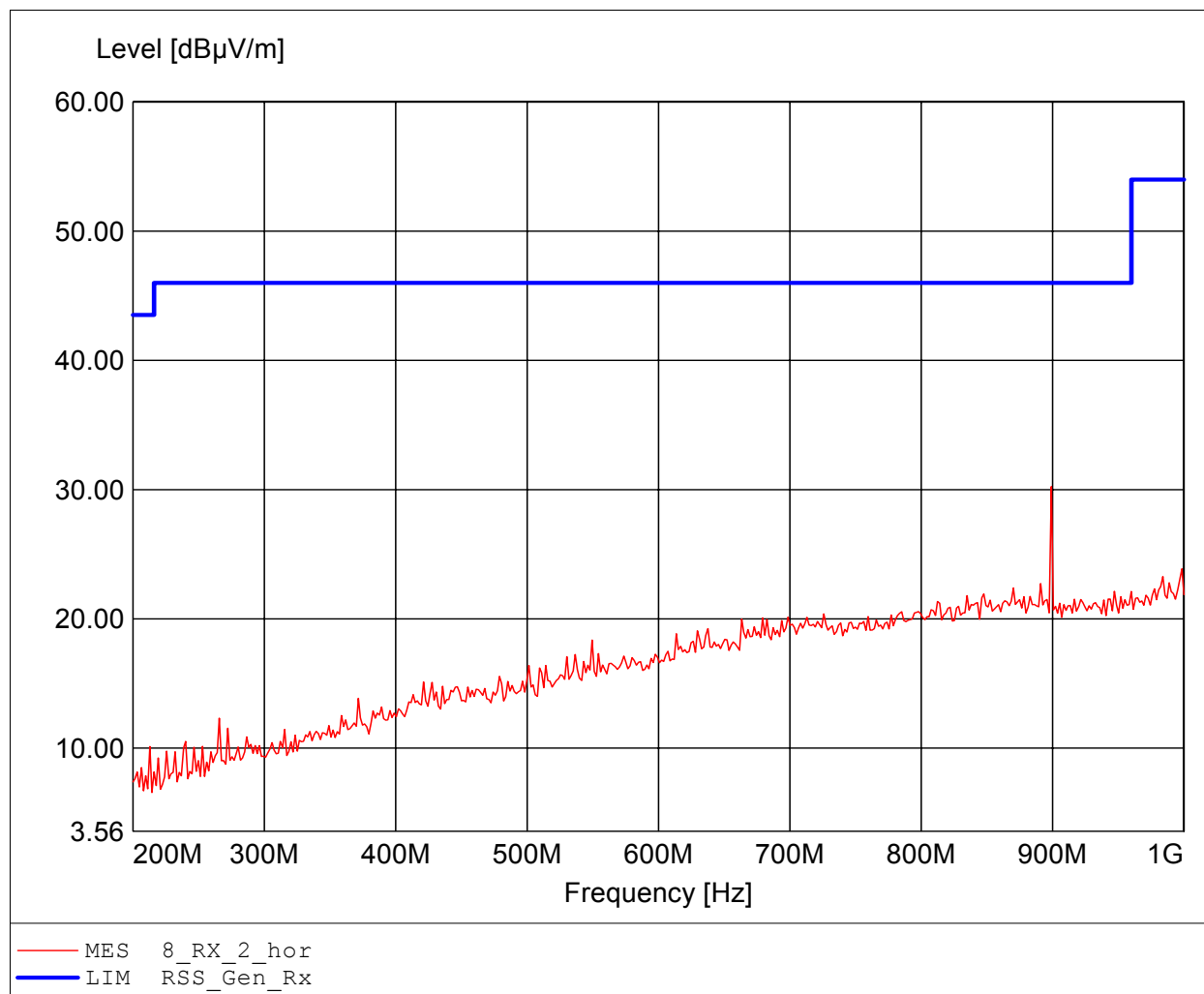

Field Strength under normal conditions

Standards Industry Canada, RSS-GEN


Approval Holder: LifeScan Scotland Ltd. / G0M-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 8 / Freq: 911MHz / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 8
Comment 1: Dist.: 3m, Ant.: HL025, ampl.
Comment 2: Freq: 3.964GHz Emax: 43.83dB μ V/m RBW: 1 MHz

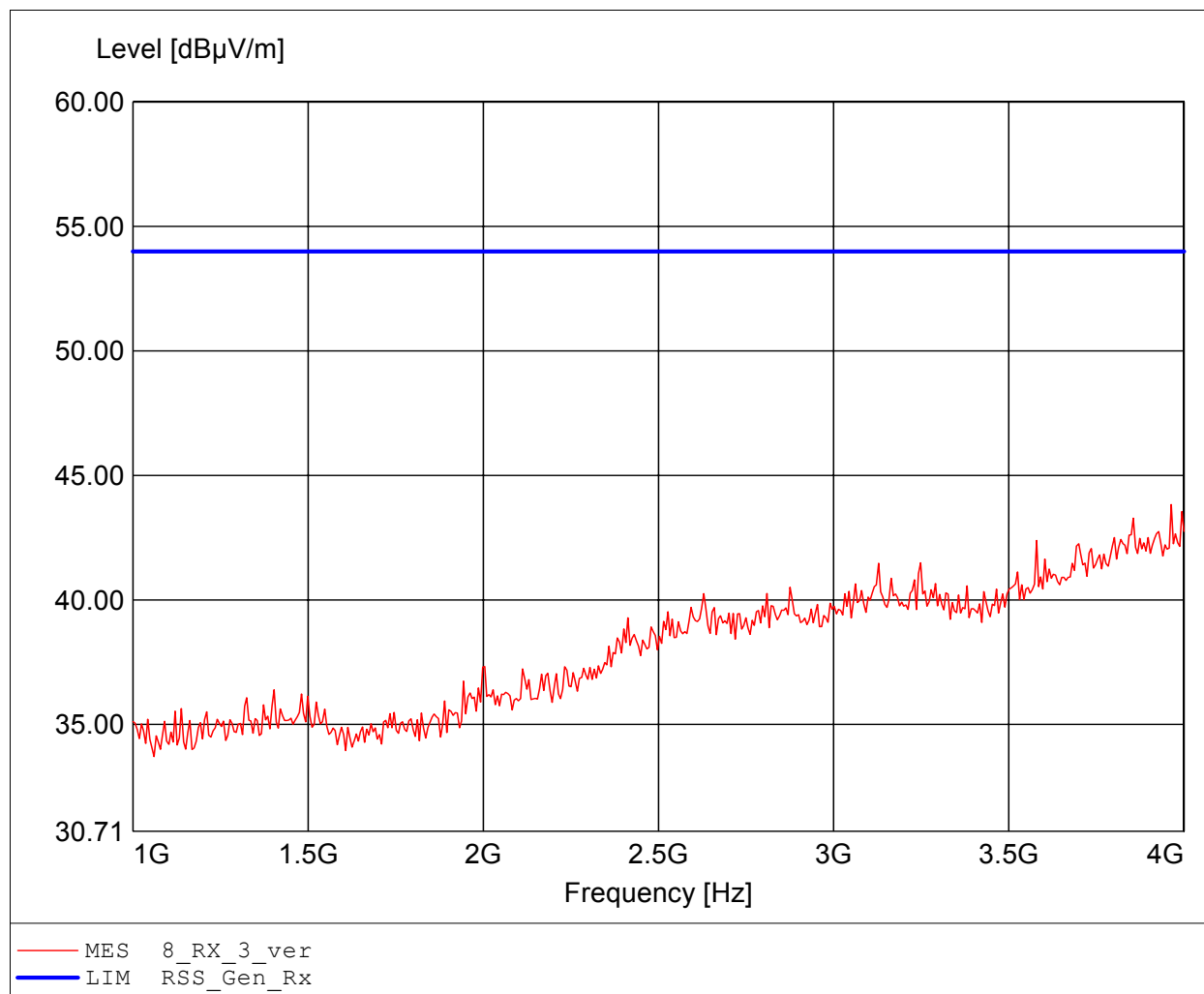

Field Strength under normal conditions

Standards Industry Canada, RSS-GEN

Approval Holder: LifeScan Scotland Ltd. / G0M-1104-1064
EUT: Blood glucose meter / Titan
Config: CH: 8 / Freq: 911MHz / vertical
Test Site / Operator: Eurofins Product Service GmbH / Mr. Pudell
Test Condition: Tnom.: 24°C / Unom: 2x 1.5 V DC (Battery)
Test Specification: Freq. / CH: 8
Comment 1: Dist.: 3m, Ant.: HL025, ampl.
Comment 2: Freq:3.934GHz Emax:42.94dBµV/m RBW: 1 MHz

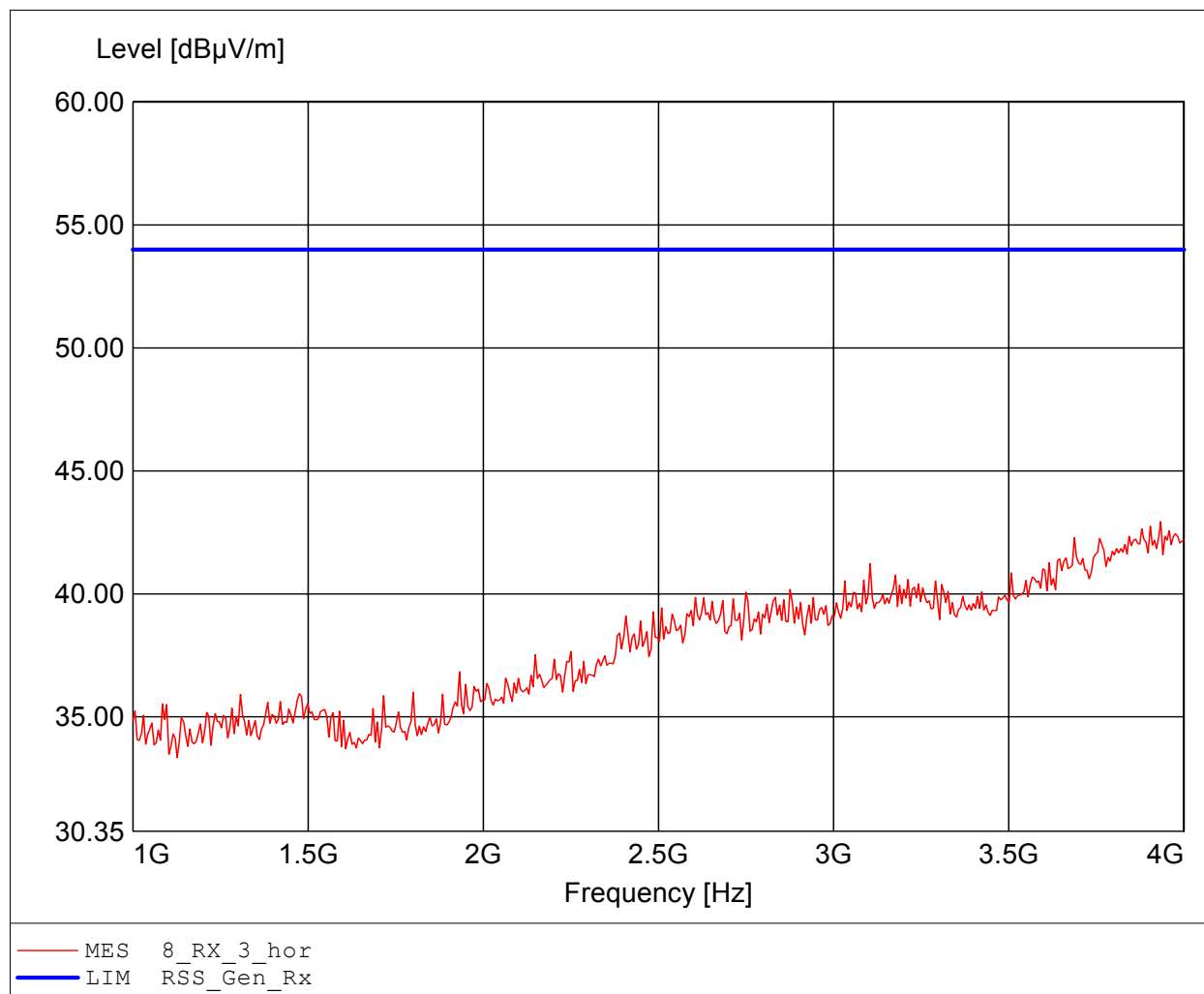


Annex F AC Power line Conducted Emissions

EMI voltage test in the ac-mains according to FCC Part 15b

Order number: G0M-1104-1064

Manufacturer: LifeScan California
 EUT Name: Blood glucose meter (Titan)
 Model: One Touch Ping Verio
 Test Site: Eurofins Product Service GmbH
 Operator: Mr. Pflug
 Test Conditions: Tnom: 23°C, Unom: 120VAC(AC/DC-adapter)
 LISN: ESH2-Z5 N
 Mode: data-link
 Test Date: 09.05.2011

Index 5

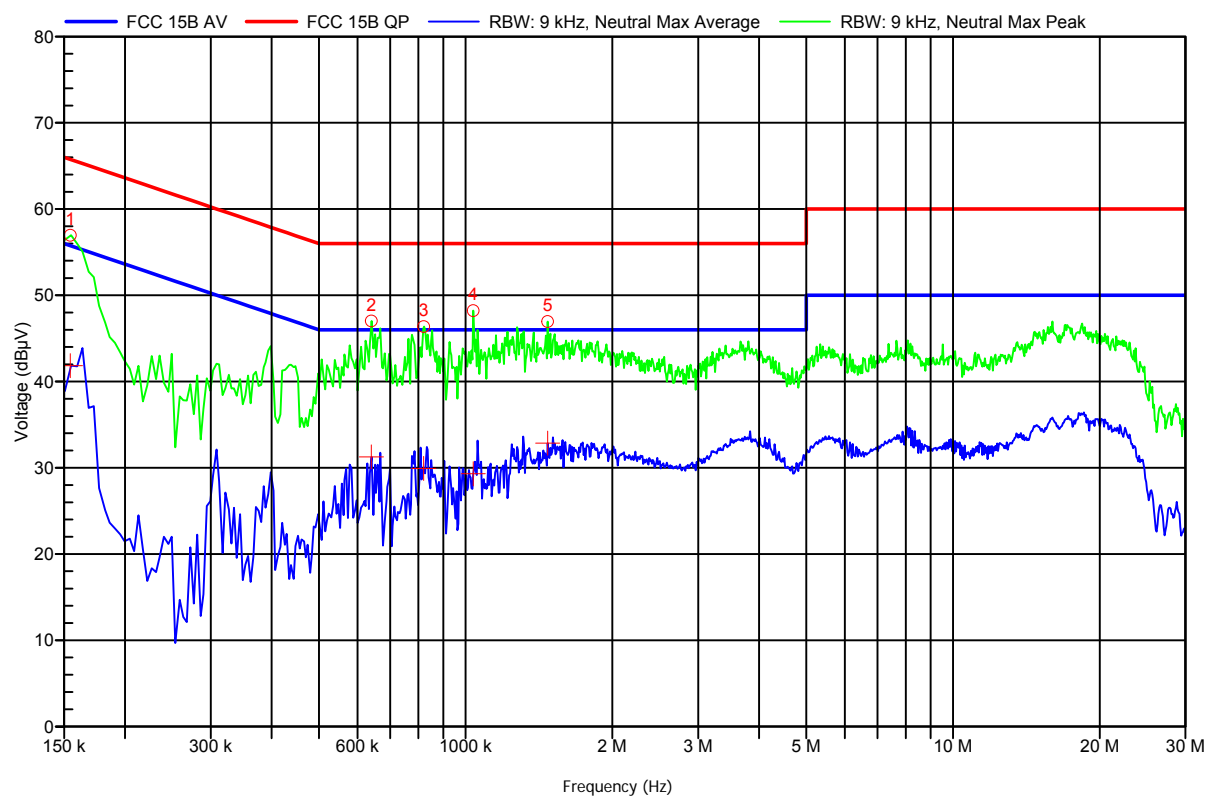

Frequency	Status
154.5 kHz	Pass
640.5 kHz	Pass
820.5 kHz	Pass
1.036 MHz	Pass
1.473 MHz	Pass

Frequency	Average	Average Limit	Average Difference	Status
154.5 kHz	41.84 dBµV	55.75 dBµV	-13.92 dB	Pass
640.5 kHz	31.26 dBµV	46 dBµV	-14.74 dB	Pass
820.5 kHz	29.93 dBµV	46 dBµV	-16.07 dB	Pass
1.036 MHz	29.3 dBµV	46 dBµV	-16.7 dB	Pass
1.473 MHz	32.86 dBµV	46 dBµV	-13.14 dB	Pass

Test Report No.: G0M-1104-1064-C-1

Eurofins Product Service GmbH
 Storkower Str. 38c, D-15526 Reichenwalde, Germany

EMI voltage test in the ac-mains according to FCC Part 15b

Order number: G0M-1104-1064

Manufacturer:	LifeScan California
EUT Name:	Blood glucose meter (Titan)
Model:	One Touch Ping Verio
Test Site:	Eurofins Product Service GmbH
Operator:	Mr. Pflug
Test Conditions:	Tnom: 23°C, Unom: 120VAC(AC/DC-adapter)
LISN:	ESH2-Z5 L
Mode:	data-link
Test Date:	09.05.2011

Index 6

